

# Keep your GC Column Alive: Tips and Tricks for Extending Column Lifetime



Ryan Birney and Alexander Ucci  
Application Engineers  
April 19, 2023



# Agenda

- Review – get to know your column
- Signs and symptoms of an unhealthy column
- How to know if your column is healthy?
- What can you do to bring your column back to life?
- Extending the lifetime of your column – sample preparation options
- Extending the lifetime of your column – backflush techniques

# The “Unboxing” of the GC Column

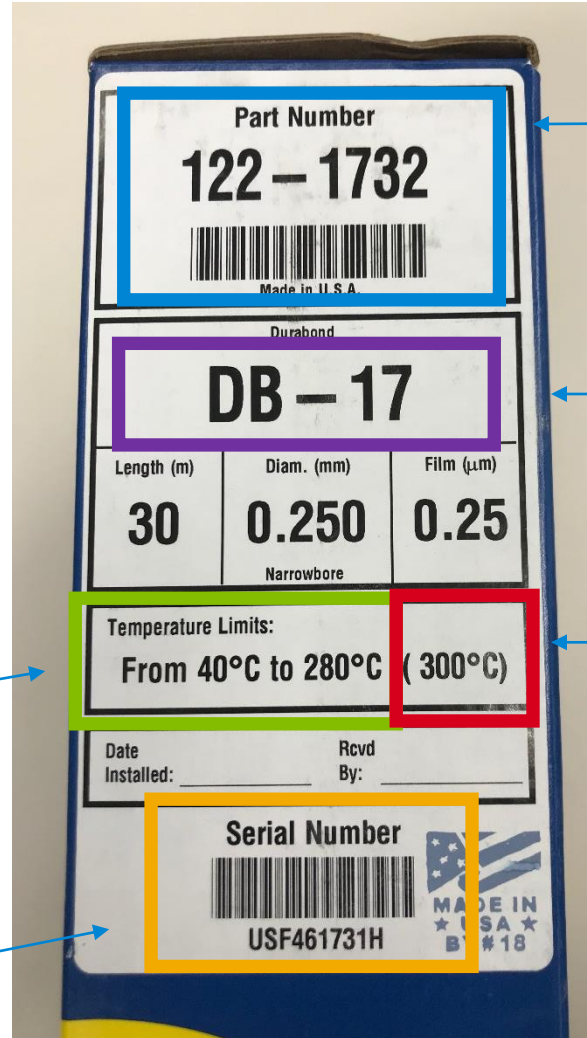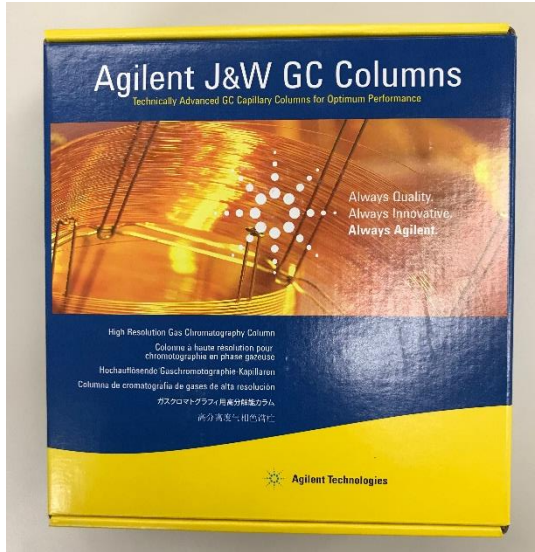

Important for identification and re-ordering


Column stationary phase

Programmed temperature limit (<10 min)

Isothermal temperature limits

Unique to each column (identification)

# What's Inside?

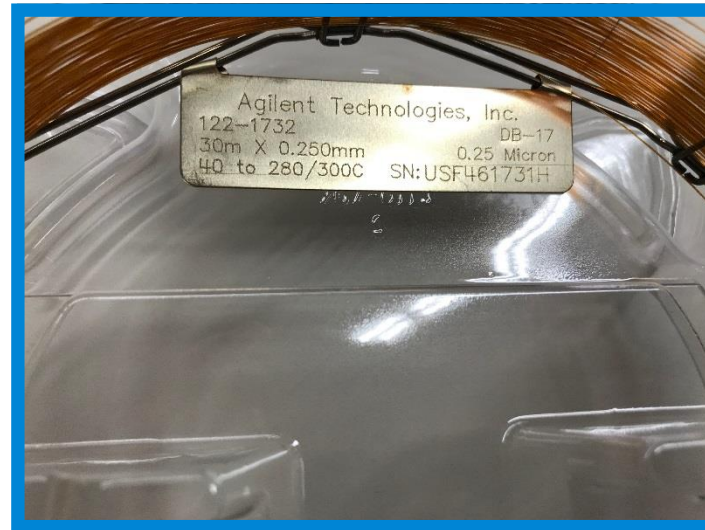

Column tag contains useful information



Column plug holds column ends together and protects against contamination. To put the column in storage, use this plug again or a piece of septa over the ends of the column.



# Column Performance Summary

Catalog: 19091S-433UI

Serial:



Stationary Phase: HP-5MS UI

Description: 30m x 0.250mm x 0.25µm

Temperature Limits: -60°C to 325°C (350°C Pgm)

## Performance Results

Theoretical Plates/Meter:

n-DECANE 3208

Retention Index:

n-PROPYLBENZENE 953.110

1-HEPTANOL 967.660

Resolution:

1-OCTENE, n-OCTANE 2.97

## Compound Identification

Compound Identification	Retent. Time	Part. Ratio	1/2-Width
1. PROPIONIC ACID	1.543	0.30	0.027
2. 1-OCTENE	2.203	0.86	0.015
3. n-OCTANE	2.282	0.92	0.016
4. 1,3-PROPANEDIOL	2.552	1.15	0.020
5. 4-METHYLPYRIDINE	3.051	1.57	0.021
6. n-NONANE	3.738	2.15	0.027
7. TRIMETHYLPHOSPHATE	4.482	2.78	0.033
8. n-PROPYLBENZENE	5.193	3.38	0.038
9. 1-HEPTANOL	5.682	3.79	0.041
10. 3-OCTANONE	6.368	4.37	0.047
11. n-DECANE	6.940	4.85	0.053

## Test Conditions


Inlet: Split (250°C) Detector: FID (325°C)

Carrier Gas: Hydrogen Flow: 42.1 cm/sec (1.2 ml/min)

Holdup Compound: Pentane (1.187-min)

Temperature Program: Isothermal at 65°C

# Chromatographic Performance

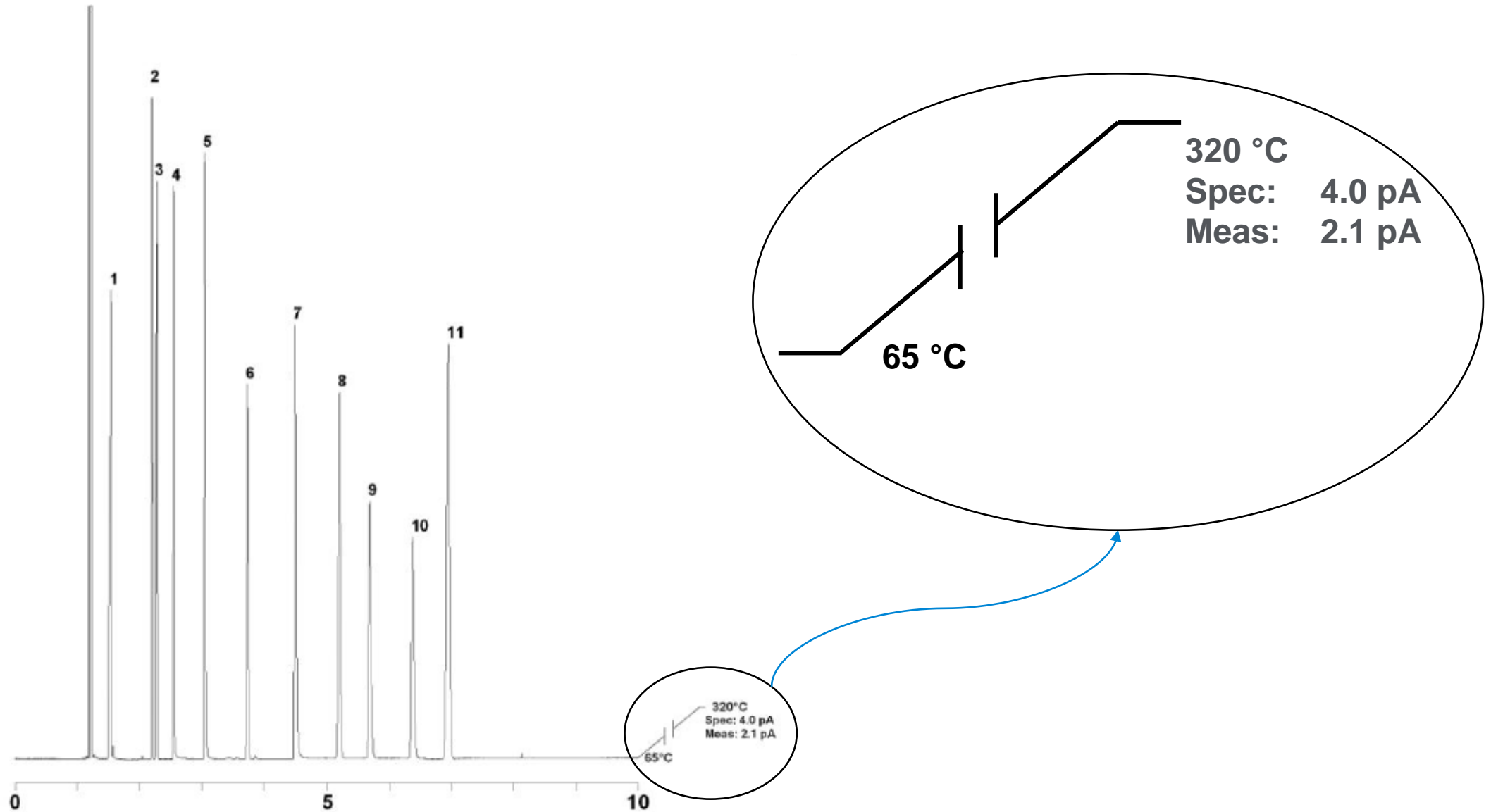

# Test Mixture Components


Compounds  
Hydrocarbons

Purpose  
Efficiency  
Retention

FAMEs, PAHs  
Alcohols  
Acids  
Bases

Retention  
Activity  
Acidic character  
Basic character

# Generating a Bleed Profile

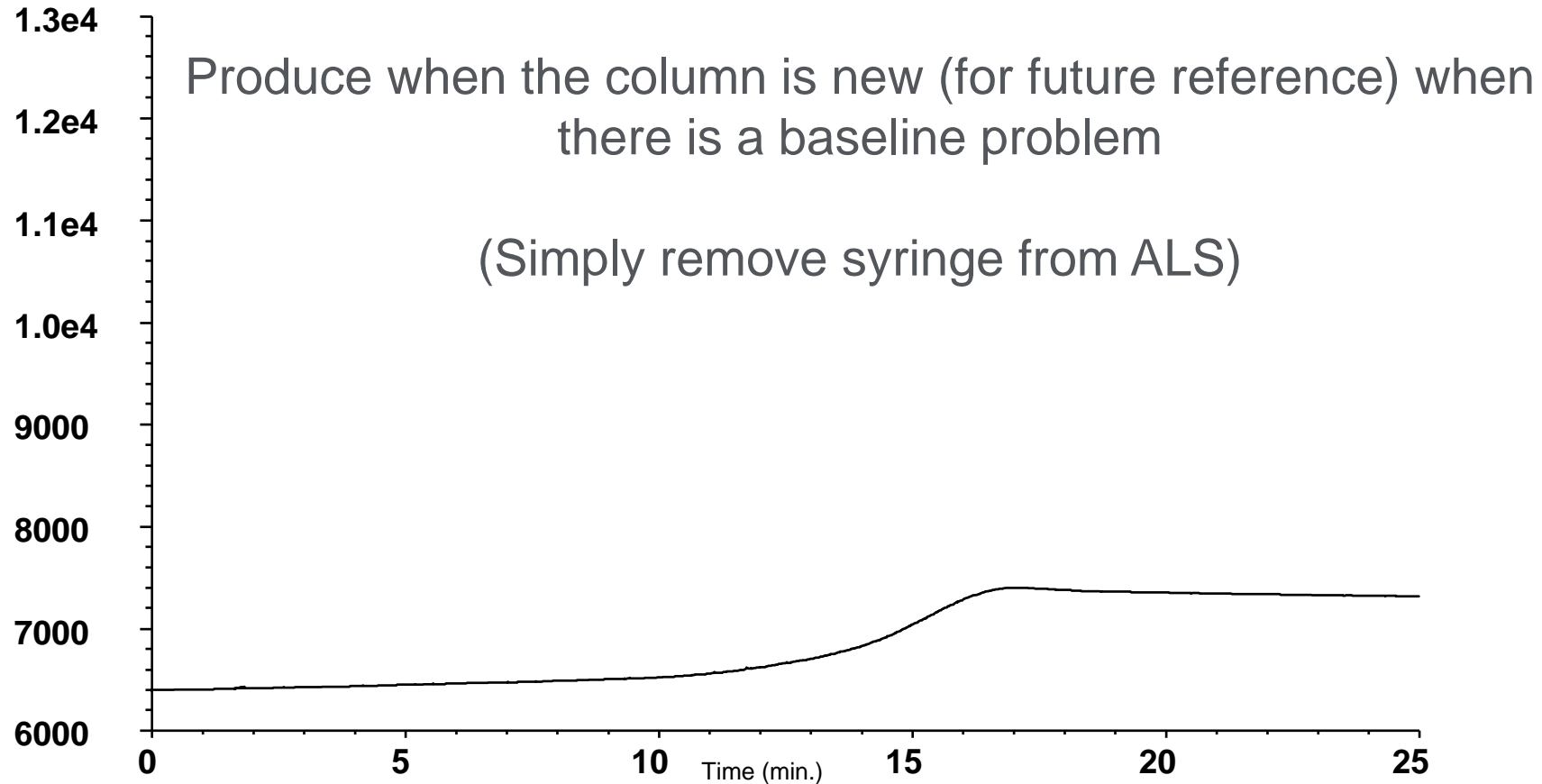

\*Agilent J&W DB-1 30 m x 0.32 mm id, 0.25  $\mu$ m  
Temperature program // 40 °C, hold 1 min // 20 °C/min to 320 °C, hold 10 min



# Common Peak Shape Issues

- **Peak tailing** – flow path or activity
- **Bonus peaks** – in sample or backflash (carry over)
- **Split peaks** – injector problems, mixed solvent
- **No peaks** – wasn't introduced, wasn't detected
- **Response changes** – activity, injector discrimination, detector problem
- **Peak fronting** – overload or solubility mismatch, injector problems
- **Shifting retention** – leaks, column aging, contamination or damage
- **Loss of resolution** – separation decreasing, peak broadening
- **Baseline disturbances** – column bleed, contamination, electronics
- **Noisy or spiking baseline** – electronics or contaminated detector
- **Quantitation problems** – activity, injector, or detector problems


# Common Causes of Column Performance Degradation

- Physical damage to the polyimide coating
- Thermal damage
- Oxidation (O<sub>2</sub> damage)
- Chemical damage by samples
- Contamination



# Physical Damage to the Polyimide Coating

- The smaller the tubing diameter, the more flexible it is
- Avoid scratches and abrasions
- Immediate breakage does not always occur upon physical damage

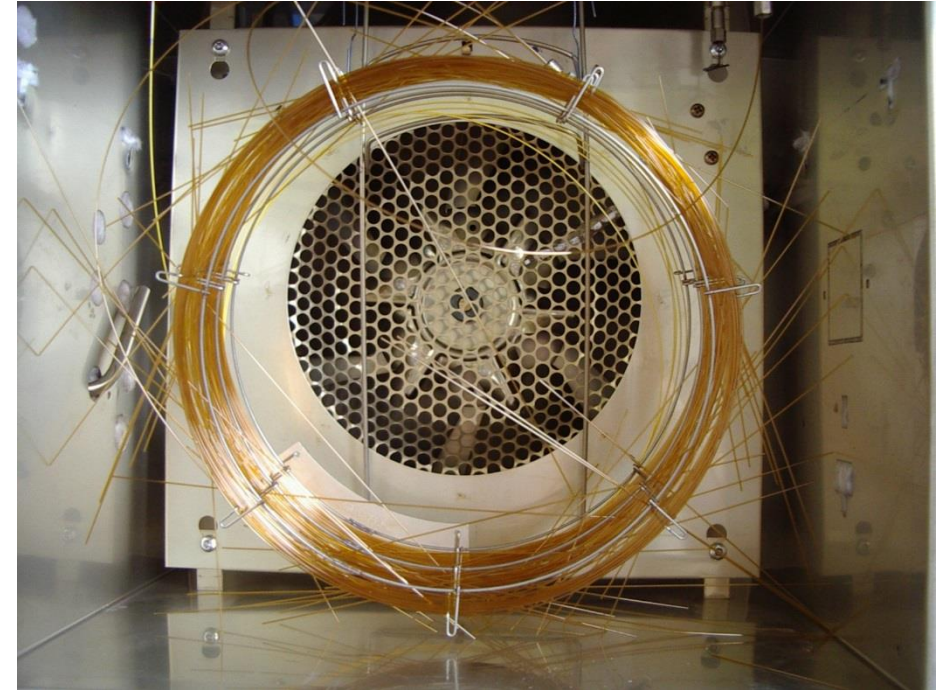

# Thermal Damage


Degradation of the stationary phase is increased at higher temperatures

- Rapid degradation of the stationary phase (breakage along the polymer backbone) caused by excessively high temperatures

Isothermal limit = indefinite time

Programmed limit = 5–10 minutes


- Temporary "column failure" below lower temperature limit
- If this happens:
  - Disconnect column from detector
  - "Bake out" overnight at isothermal limit
  - Remove 10–15 cm from column end

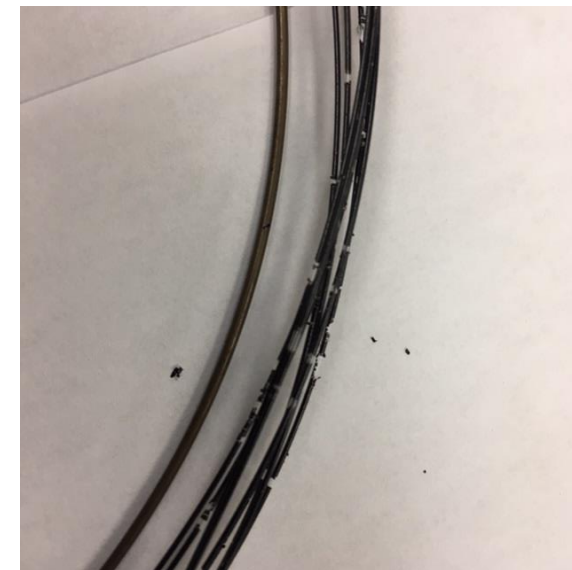



Column continuously exposed to temperatures above its temperature limit

# Oxidation (O<sub>2</sub> Damage)

Oxygen in the carrier gas rapidly degrades the stationary phase. The damage is accelerated at higher temperatures. Damage along the polymer backbone is irreversible. (Premature filament failure/excessive source maintenance.)

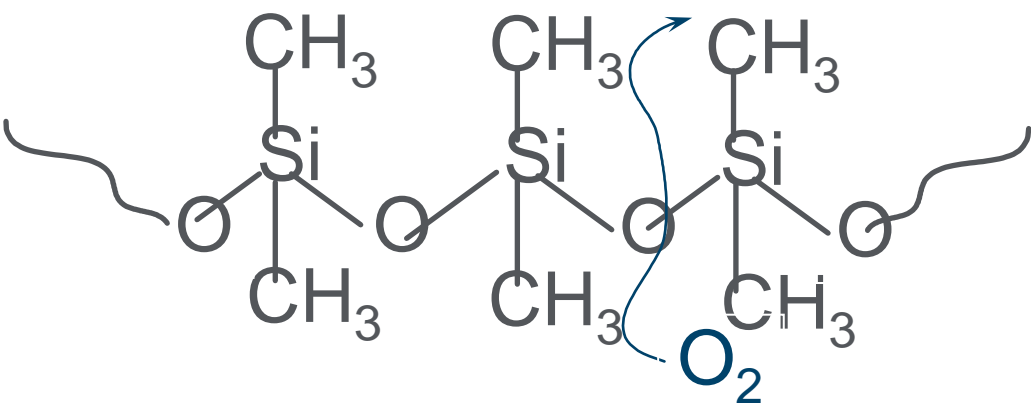


Dimethylpolysiloxane

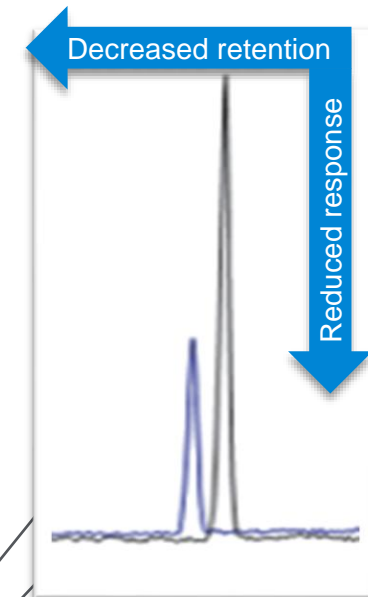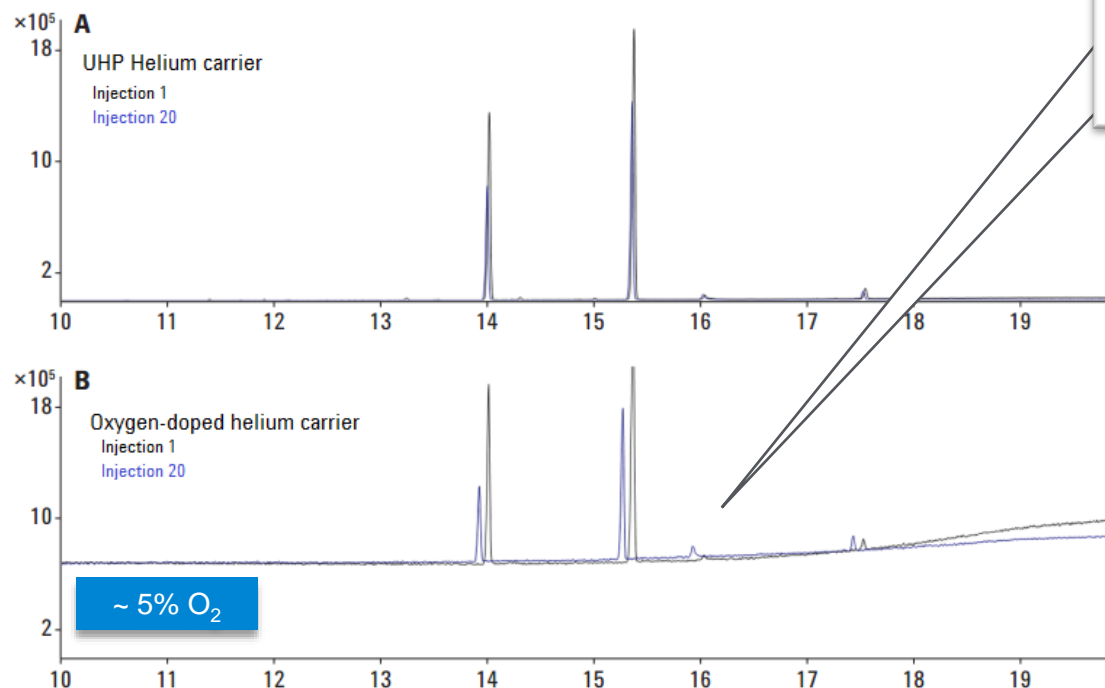


# How to Prevent Column Damage by Oxygen

- High-quality carrier gas (four 9s or greater)
- Leak free injector and carrier lines
  - Change septa
  - Maintain gas regulator fittings
- Appropriate impurity traps

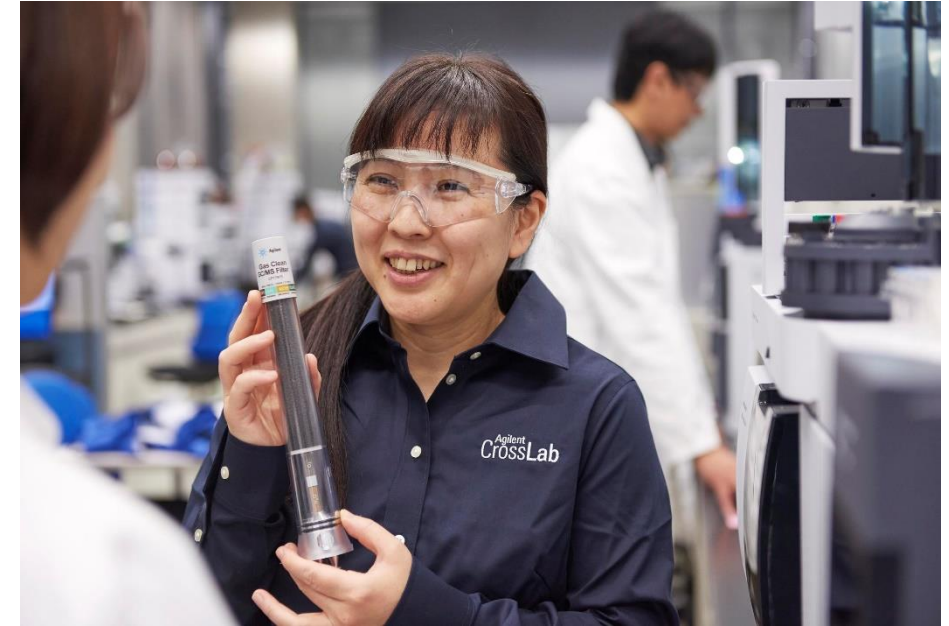

Efficient, fast, easy

# Effect of Oxygen on Peak Shape of 2-ethylhexanoic Acid

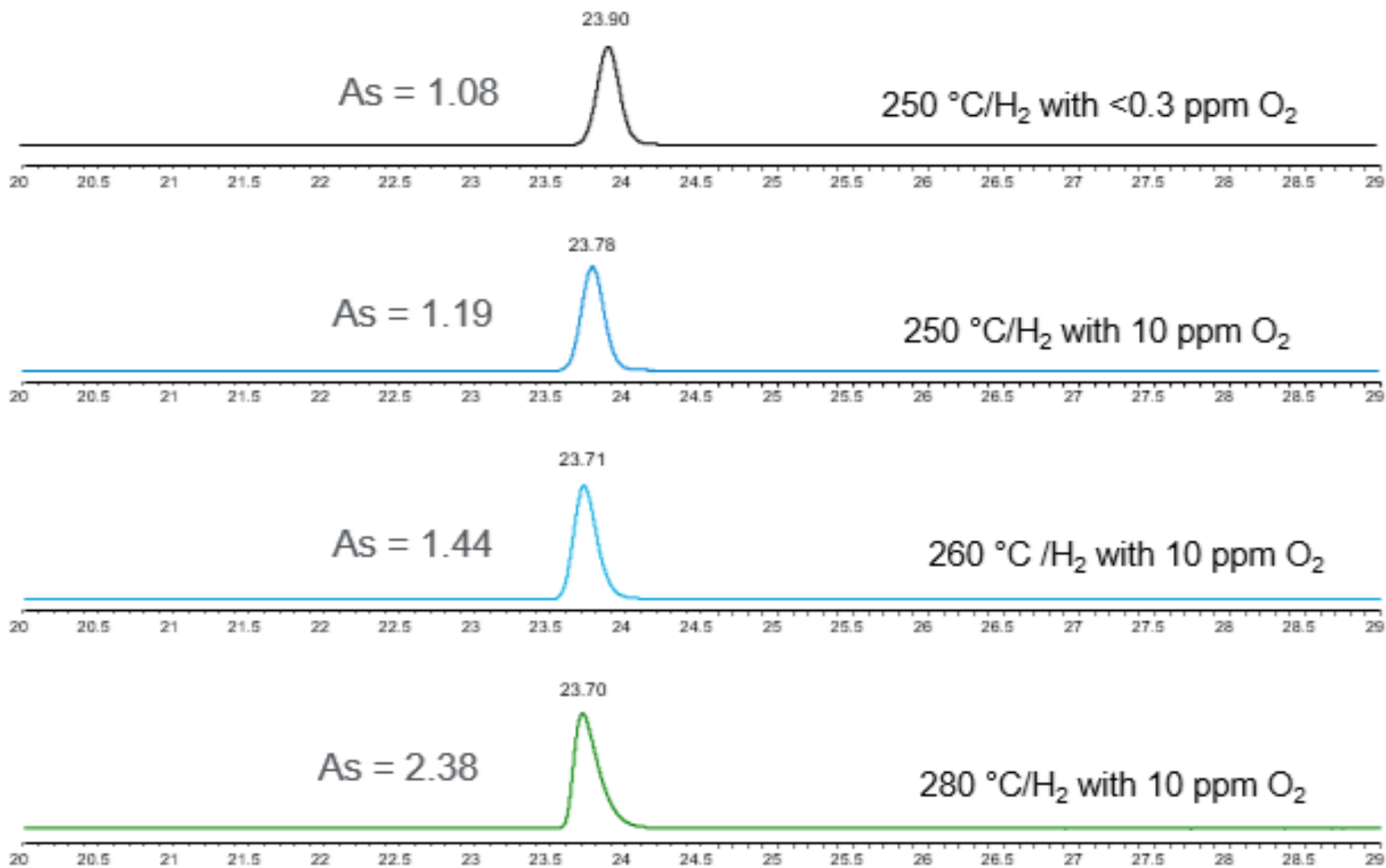

# Self Tightening Nuts: No Leaks, No Downtime, No Frustration



- Spring-driven piston continuously presses against ferrule
- Automatically retightens when ferrule shrinks
- Wing design for finger tightening
- No tools needed
- Works only with graphite/Vespel ferrules


Part Number	Description
G3440-81013	Column Nut, Collared Self-Tightening MSD
G3440-81011	Column nut, Collared Self Tightening Inlet/Detect
G3440-81012	Collar for Self Tightening Nut

<https://www.agilent.com/en/video/gc-supplies-innovation>


<https://www.agilent.com/en/video/stcn-inlet-detector>

<https://www.agilent.com/en/video/stcn-mass-spec>


Standard column nuts new fitting

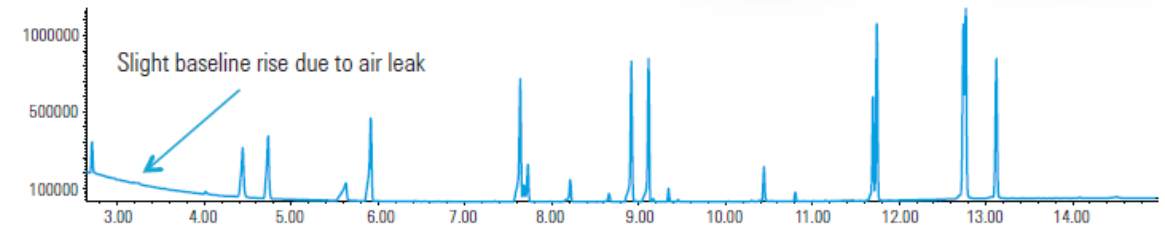

Standard column nuts after 25 injections

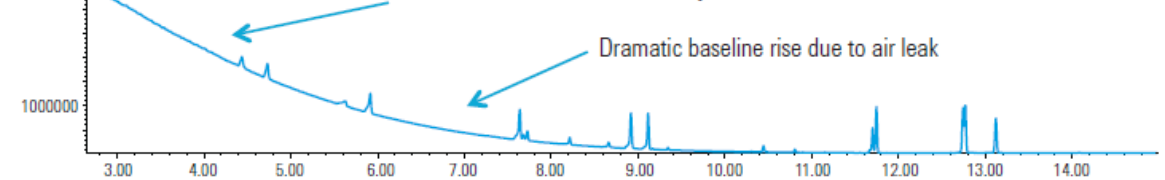

Agilent Self Tightening Column Nuts after 400 injections

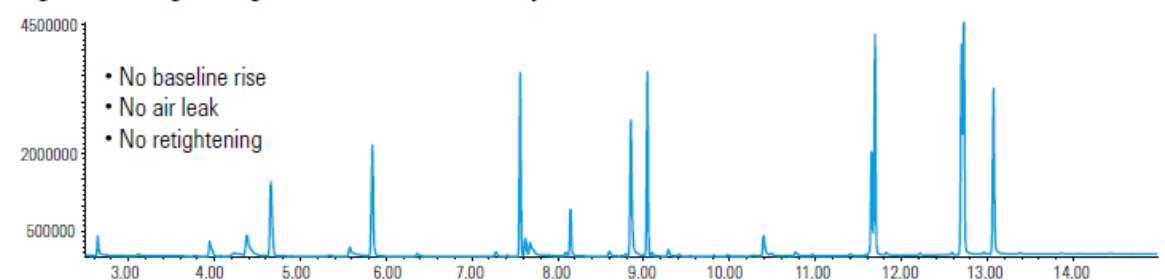

400 injections



# New Agilent Standard Winged Nut and Depth Guide

- Compatible with Agilent/HP style compact ferrules, including graphite ferrules
- Winged fastener design for easy engagement and tool free install
- Hollow-body design with low thermal mass mitigates thermal lag during temperature cycling within the GC oven
- Removable locking-collar with soft-PTFE insert to secure column placement during install without damaging the analytical column



- Easy-to-use template provides critical capillary column installation for the most popular Agilent GC-configurations
  - SSL, MMI, purge-packed inlets
  - FID, TCD, NPD detectors
  - EI MSD source
- Compatible with the Agilent Self Tightening and winged column nuts



Agilent Capillary Column Depth Guide – G3440-88000

# The Cost of Leaks

- Cost of gases
- Contamination from exposure
- Reduced consumable lifetime
- Reduced productivity from downtime
- Detector noise and elevated baselines
- Time in troubleshooting

It is critical that every customer checks for leaks. They should have the best tool for the job.

Check valves, fittings, and traps for leaks after every **piece of** maintenance, and after thermal cycling, as these can loosen some types of fittings.

# Agilent CrossLab CS (Cartridge System)

No peaks from leaks

Features:

- Exchangeable cartridge with ADM Flow Meter
- Automatic notification of probe filter replacement
- Ergonomic and robust design
- Universal 3AA or USB power
- USB connects to web interface for added functionality and firmware updates
- Easy to view OLED screen
- Kickstand

Leak detector  
cartridge

Handheld



ADM Flow Meter  
cartridge

# More Information

- **Agilent.com CrossLab CS Leak Detector**

[www.agilent.com/chem/gas-leak-detector](http://www.agilent.com/chem/gas-leak-detector)

- **Agilent.com – ADM Flow Meter**

<https://www.agilent.com/en/product/gas-purification-gas-management/gas-management/adm-flow-meter>

- **Installation manual**

*Agilent CrossLab CS  
Electronic Leak Detector manual*

Part number: G6693-90000

The installation manual is available on Agilent.com.

- **Innovation minute video**

<https://www.agilent.com/en/video/crosslab-cs-innovation-minute>

- **Technical overview**

*Agilent CrossLab Cartridge System  
(CS) Electronic Leak Detector*

Publication number: 5994-4262EN

The technical overview is available on Agilent.com.

- **Brochure**

*GC Troubleshooting in  
the Palm of Your Hand*

Publication number: 5994-3607EN

The brochure is available on Agilent.com.

- **Flyer**

*Is a Leak Causing Your  
Inaccurate Results?*

Publication number: 5994-4202EN

The flyer is available on Agilent.com.

# Ordering Guide

## One year warranty

- G6693A – CrossLab CS Electronic Leak Detector
- G6694A – Electronic Leak Detector Cartridge
- G6699A - CrossLab CS Bundle: ADM Flow Meter and Electronic Leak Detector
  - The bundle will include one handheld, two cartridges, and a **free** carrying case.
- G6694-60005 – Replacement probe filter
- G6691-40500 – Carrying case



## Existing products:

- G6691A – CrossLab CS ADM Flow Meter
- G6692A – ADM Flow Meter cartridge\*
- Note that the ADM Flow Meter cartridge is ordered annually for calibration. The Electronic Leak Detector does not need to be recalibrated

# Chemical Damage

Bonded and crosslinked columns have excellent chemical resistance except for inorganic acids and bases.

HCl NH<sub>3</sub> KOH NaOH

H<sub>2</sub>SO<sub>4</sub> H<sub>3</sub>PO<sub>4</sub> HF

Chemical damage will be evident by excessive bleed, lack of inertness, or loss of resolution/retention.

# Chemical Damage

## What to do if it happens


- Remove 0.5 to 1 m from the front of the column
- Severe cases may require removal of up to 5 m



# What is Normal Column Bleed?

Normal background signal generated by the elution of normal degradation products of the column stationary phase. Column bleed is influenced by:


- Phase type
- Temperature
- Column dimensions

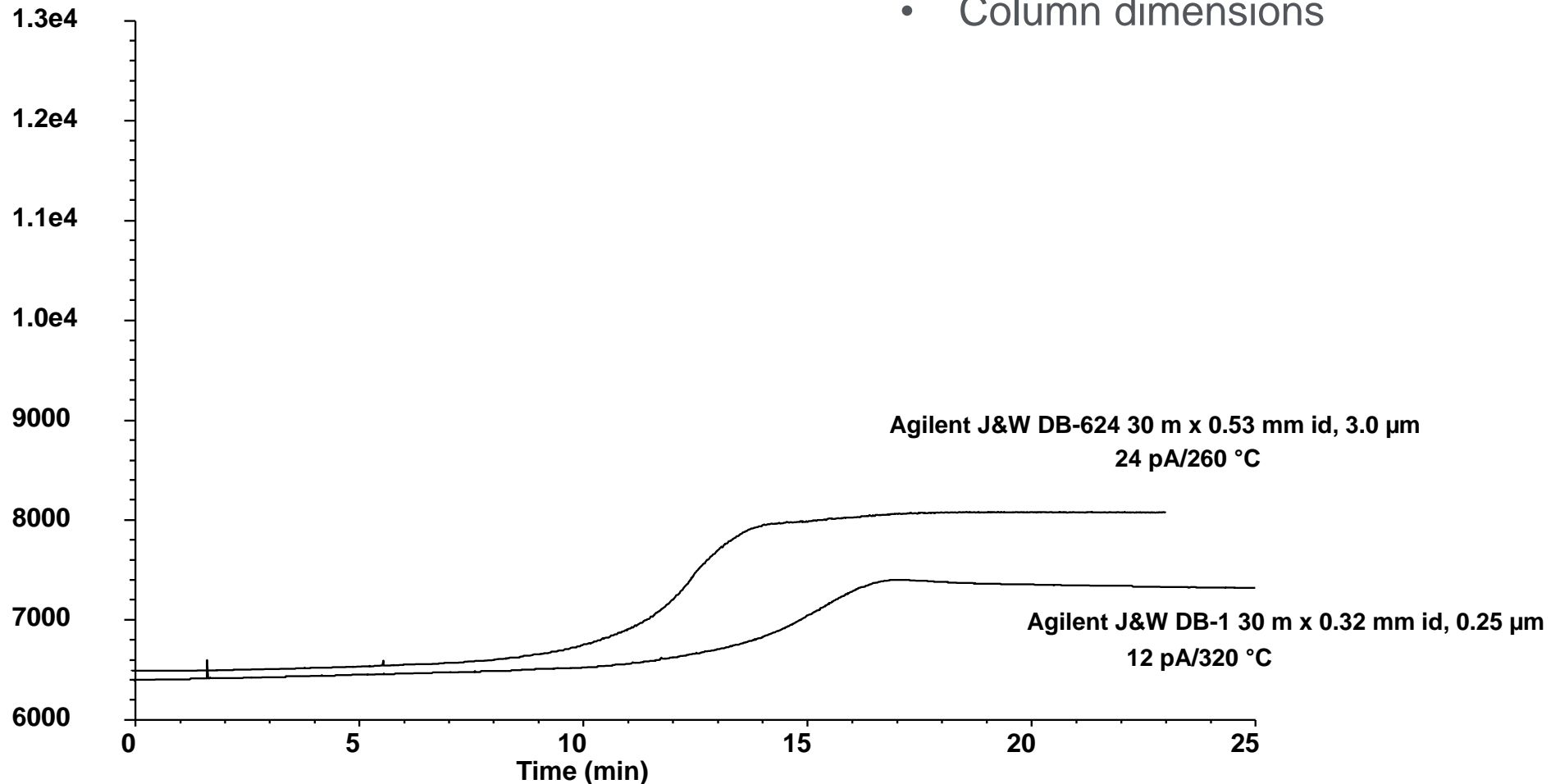



# Mass Spectrum of Phenylmethylpolysiloxane Column Bleed

Normal background (HP-5ms UI)

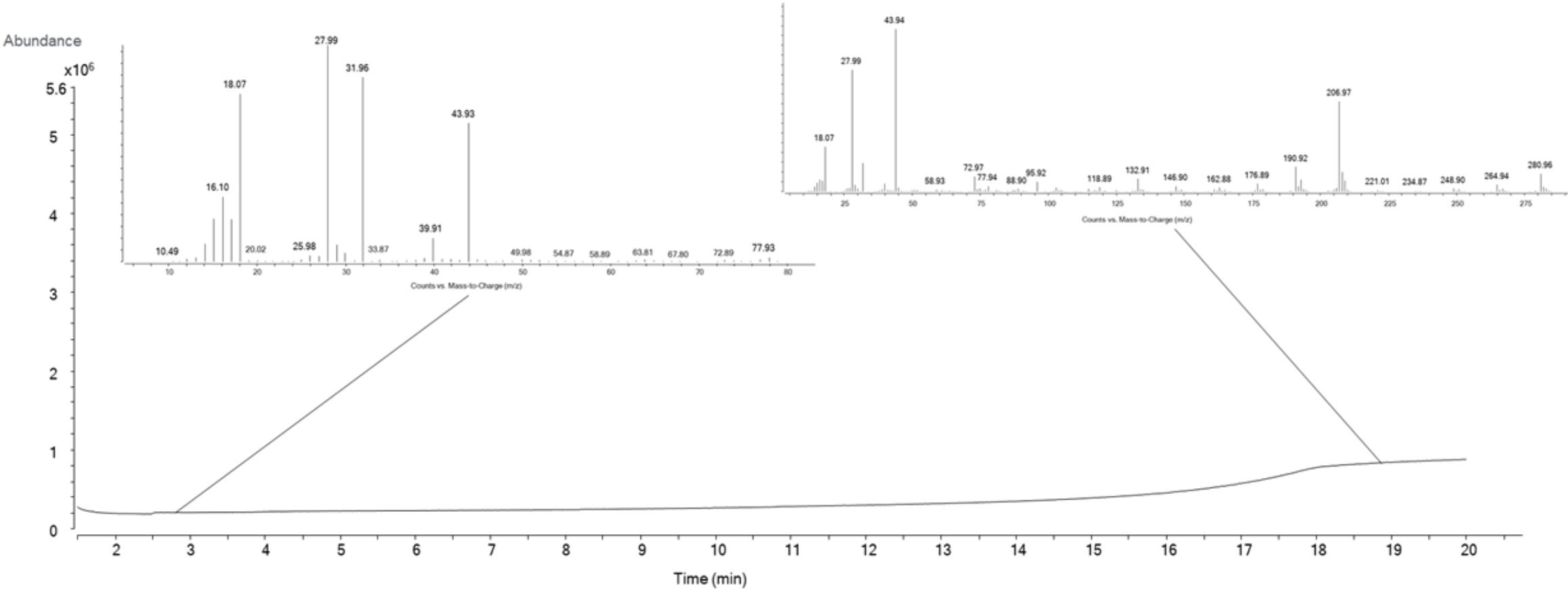

# What is a Bleed Problem?

An abnormal elevated baseline at high temperature

It is **not**




- A high baseline at low temperature
- Wandering or drifting baseline at any temperature
- Discrete peaks

# Column Contamination and Symptoms

- Fouling of GC and column by contaminants
- Mimics nearly every chromatographic problem
- Poor peak shape
- Loss of separation (resolution)
- Changes in retention
- Reduced peak size
- Baseline disturbances (semivolatiles only)



# Example of Column Contamination

Agilent J&W DB-624 QC test mix\*  
After 75 injections of oily sample

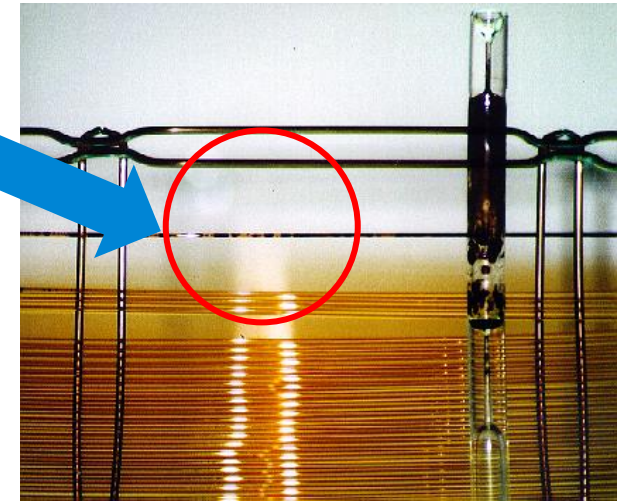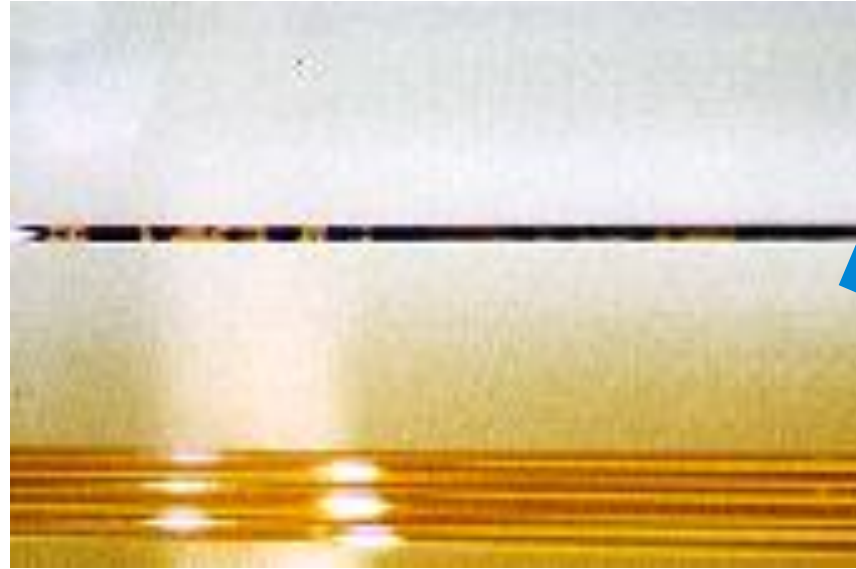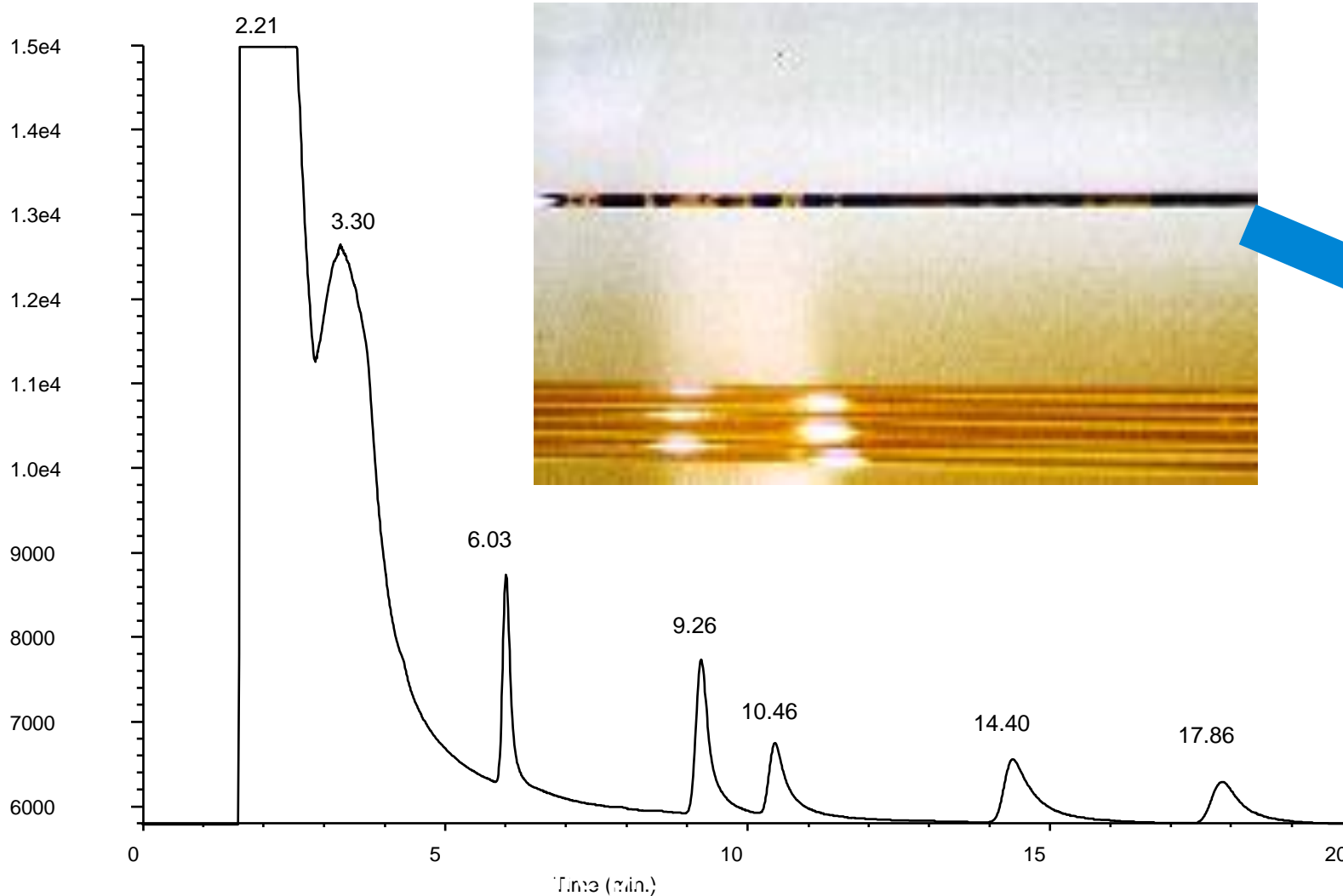


# Sources of Contamination

- Septum and ferrule particles
- Gas and trap impurities
- Unknown sources (vials, syringes, and so on)

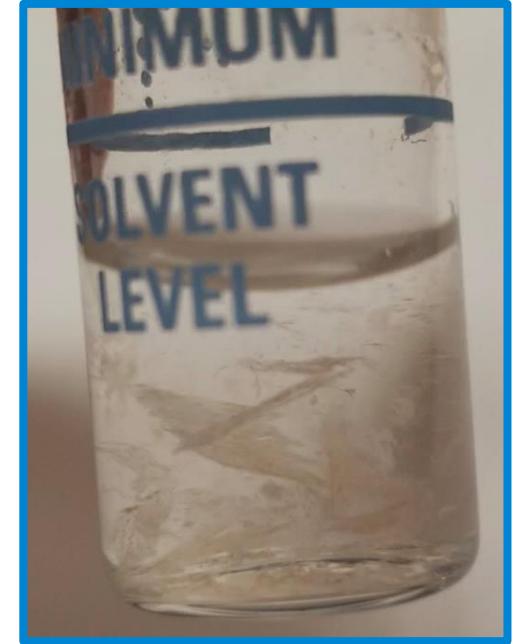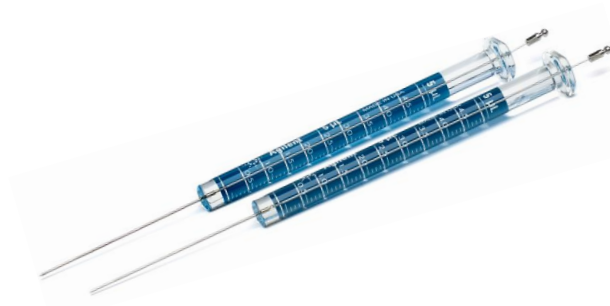

Sample vial septum bleed profile:

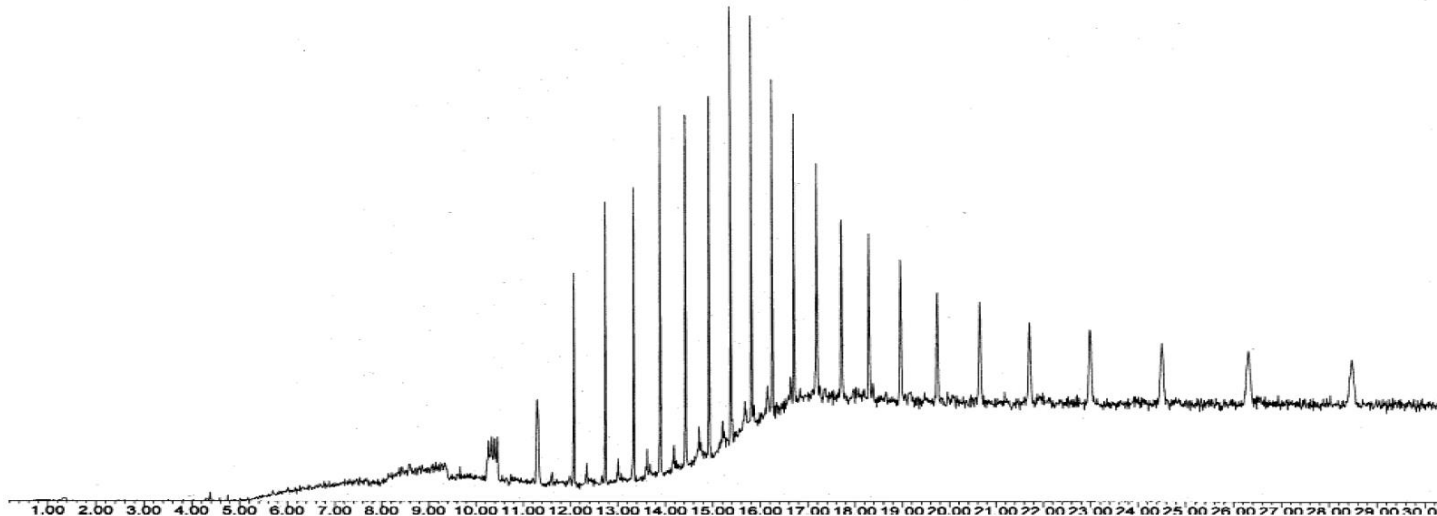

Contaminated wash solvent

# Typical Samples That Contain a Large Amount of Residues

Biological  
Soils  
Foods  
Wastewater  
Sludges

All samples contain residues (even standards)



# Types of Residues

## Nonvolatile residues

- Any portion of the sample that does not elute from the column or remains in the injector.

## Semivolatile residues

- Any portion of the sample that elutes from the column after the current chromatographic run.

# Semivolatile Contamination

## What to do if it happens

- “Bake out” the column
  - Limit to 1 to 2 hours
  - Longer times may polymerize some contamination and reduce column life
- Solvent rinse the column



# Nonvolatile Contamination

## What to do if it happens

- Do not “bake out” the column
- Front end maintenance
  - Clean or change the injector liner
  - Clean the injector
  - Cut off 0.5 to 1 m of the front of the column
- Turn the column around
- Cut the column in half

# Methods to Minimize Nonvolatile Residue Problems

- Sample cleanup
- Packed injection port liners
- Guard columns
- Backflush techniques



# Dilute and Shoot

## Advantages

- Fast and easy
- High throughput



GC inlet liner



GC inlet seal

## Limitations

- Interferences are not removed
- Analyte concentration is reduced
- Instrument and column contamination
- Matrix interferences – ion suppression or poor peak shapes

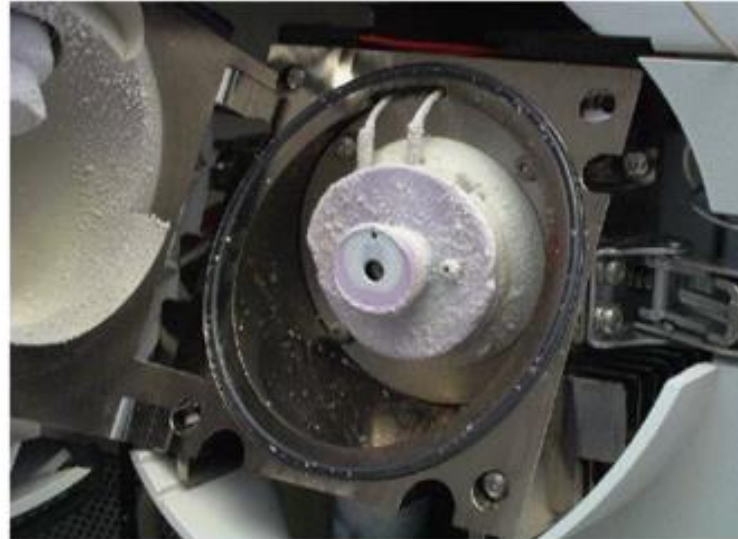


Image of salt buildup on an ESI-LC/MS inlet from unremoved salts.

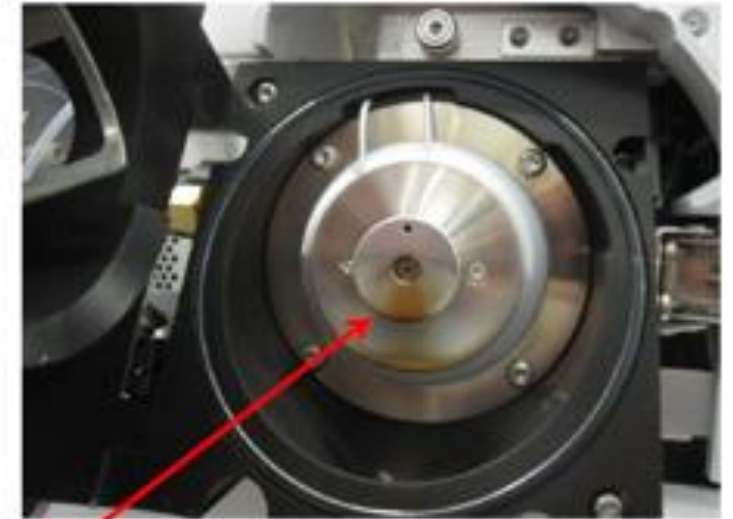


Image shows the build up on the ESI-MS inlet after 3000X urine dilute and shoot injections.

# Physical Filtration


## Captiva premium syringe filters

- Certified to be free of UV-detectable extractables on HPLC. PES and glass fiber also certified for LC/MS.
- Color-coded boxes for easy identification
- Comprehensive portfolio to meet all customers' needs

Premium Syringe Filters						
Membrane	Diameter/Pore Size					
	4 mm		15 mm		25 mm	(28 mm)
	0.2 µm	0.45 µm	0.2 µm	0.45 µm	0.2 µm	0.45 µm
PTFE	◆	◆	◆	◆	◆	◆
Nylon			◆	◆	◆	◆
PES	◆	◆	◆	◆	◆	◆
Regenerated cellulose	◆	◆	◆	◆	◆	◆
Cellulose acetate					◆	◆
Glass microfiber			◆		◆	
Depth filters: glass/PTFE			◆	◆	◆	◆
Depth filters: glass/nylon			◆	◆	◆	◆



# Physical Filtration – Captiva Filter Vials

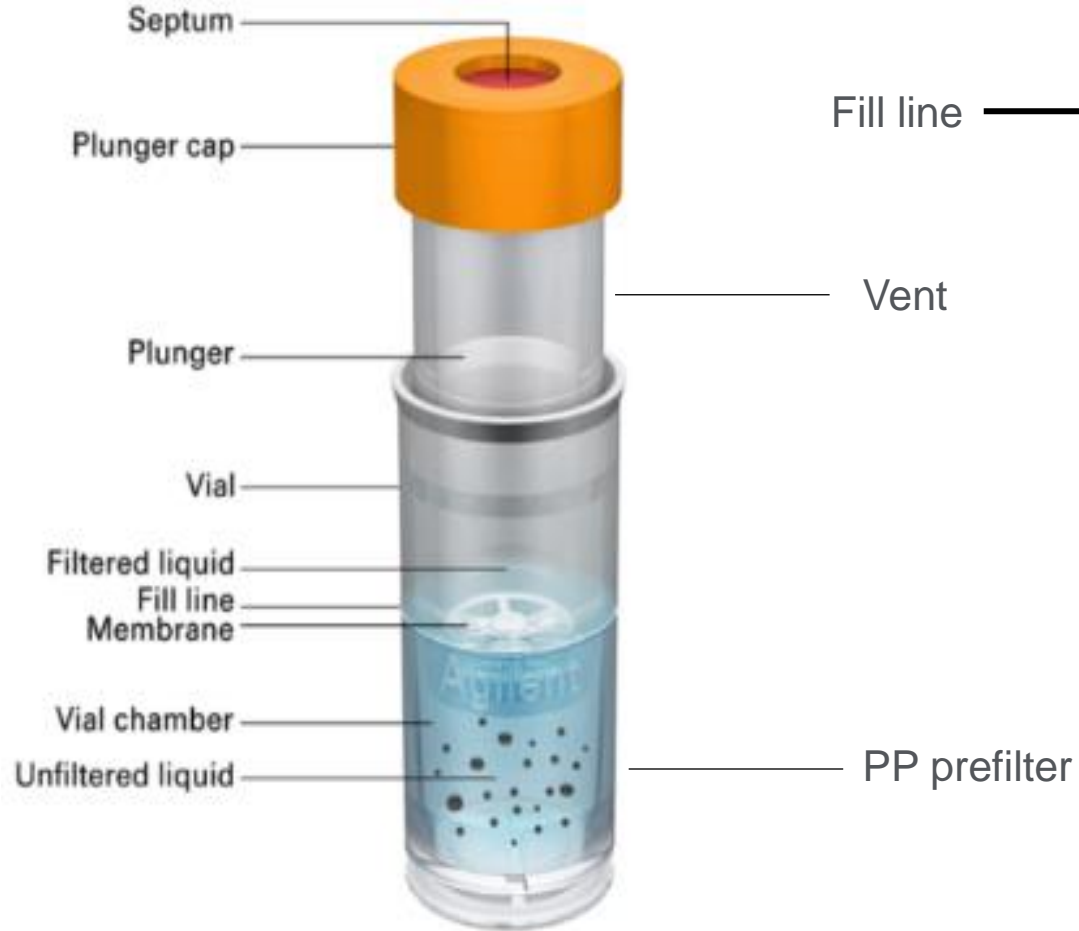

Part Number	Description
5191-5933	PTFE filter vial, 0.45 $\mu\text{m}$ , 100/pk
5191-5934	PTFE filter vial, 0.20 $\mu\text{m}$ , 100/pk
5191-5935	Nylon filter vial, 0.45 $\mu\text{m}$ , 100/pk
5191-5936	Nylon filter vial, 0.20 $\mu\text{m}$ , 100/pk
5191-5939	RC filter vial, 0.45 $\mu\text{m}$ , 100/pk
5191-5940	RC filter vial, 0.20 $\mu\text{m}$ , 100/pk
5191-5941	PES filter vial, 0.45 $\mu\text{m}$ , 100/pk
5191-5942	PES filter vial, 0.20 $\mu\text{m}$ , 100/pk
5191-5943	Vial closure tool

Preslit options available as well

Agilent.com/chem/filtervials  
Filter vials user guide: 5994-0814EN

# Chemical Filtration

## Captiva EMR-Lipid


- One of the newest Agilent sample cleanup products with a 2-in-1 benefit of removing lipids and fats selectively and efficiently.
- It reduces ion suppression, increases analyte sensitivity, improves peak shape, and extends the lifetime of your analytical column.
- Simple pass-through format, 96-well plate, 1 mL, 3 mL, and 6 mL cartridges
- Solvent-retention frit in 1 mL cartridge/96-well plate for in-well protein precipitation
- Unique chemistry and filtration ensures protein and lipid removal
- Depth filtration design allows for smooth elution
- Received the Analytical Scientist Innovation Award (TASIA) of 2017



# Captiva EMR-Lipid Sorbent Technology

EMR-Lipid sorbent technology effectively traps lipids through two mechanisms:


- **Size exclusion** – Unbranched hydrocarbon chains (lipids) enter the sorbent; bulky analytes do not
- **Sorbent chemistry** – Lipid chains that enter the sorbent are trapped by hydrophobic interactions

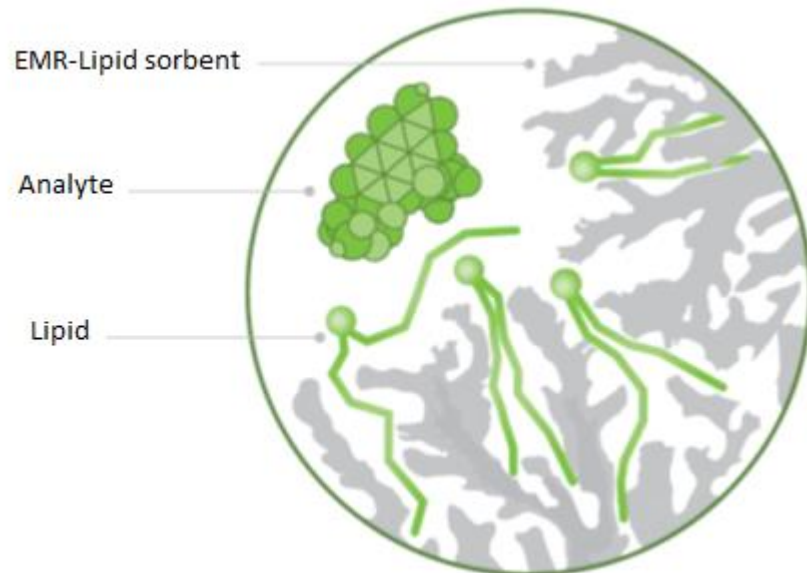




# Captiva EMR-Lipid

## Selective removal of lipids

Removes lipids

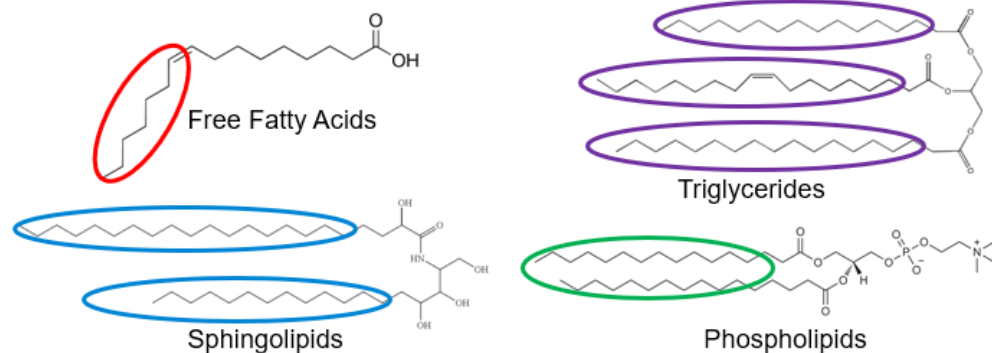




Does not remove target analytes

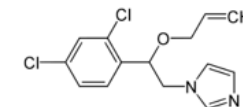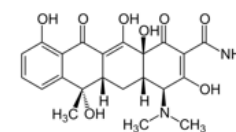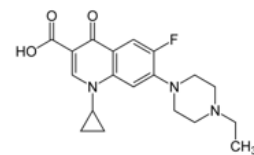

Fluoroquinolones

Tetracyclines

Imidazole pesticides

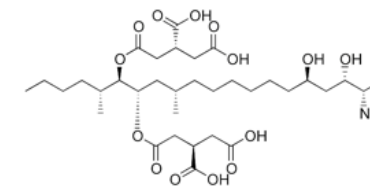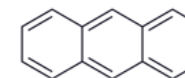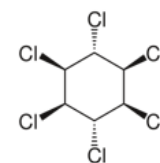

Organochlorine Pesticides

PAHs


Fumonisin B2



# Captiva EMR-Lipid


General protocol for food and food products using 3 mL and 6 mL cartridges

## Operating instructions

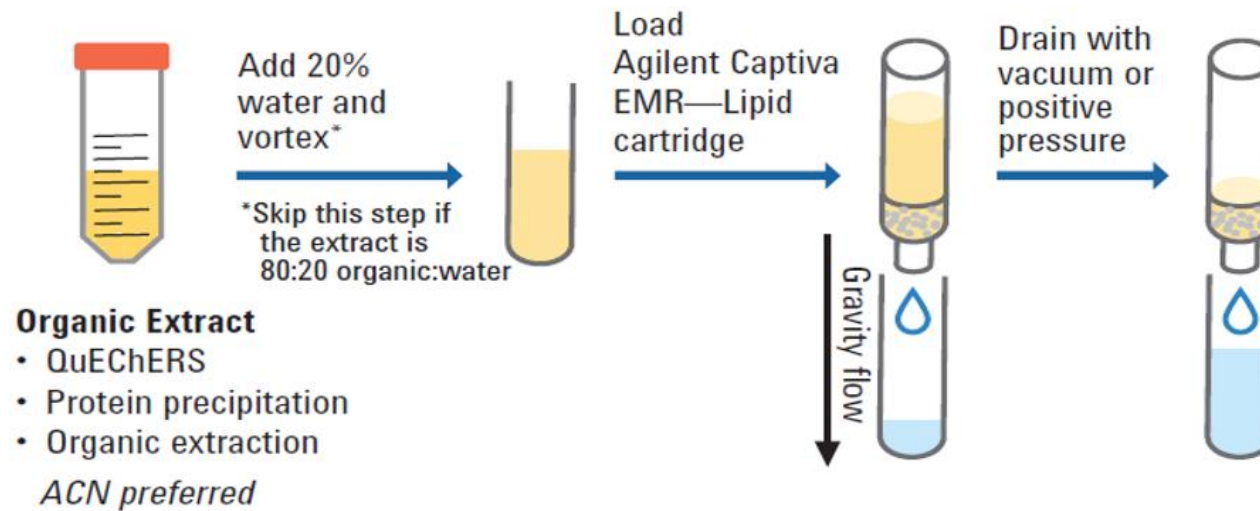


[Captiva EMR-Lipid method guide for 3 mL and 6 mL cartridges](#)

# A Cleanup Step Improves Analytes, S/N Ratio, and Integration Accuracy on GC/MS(/MS)

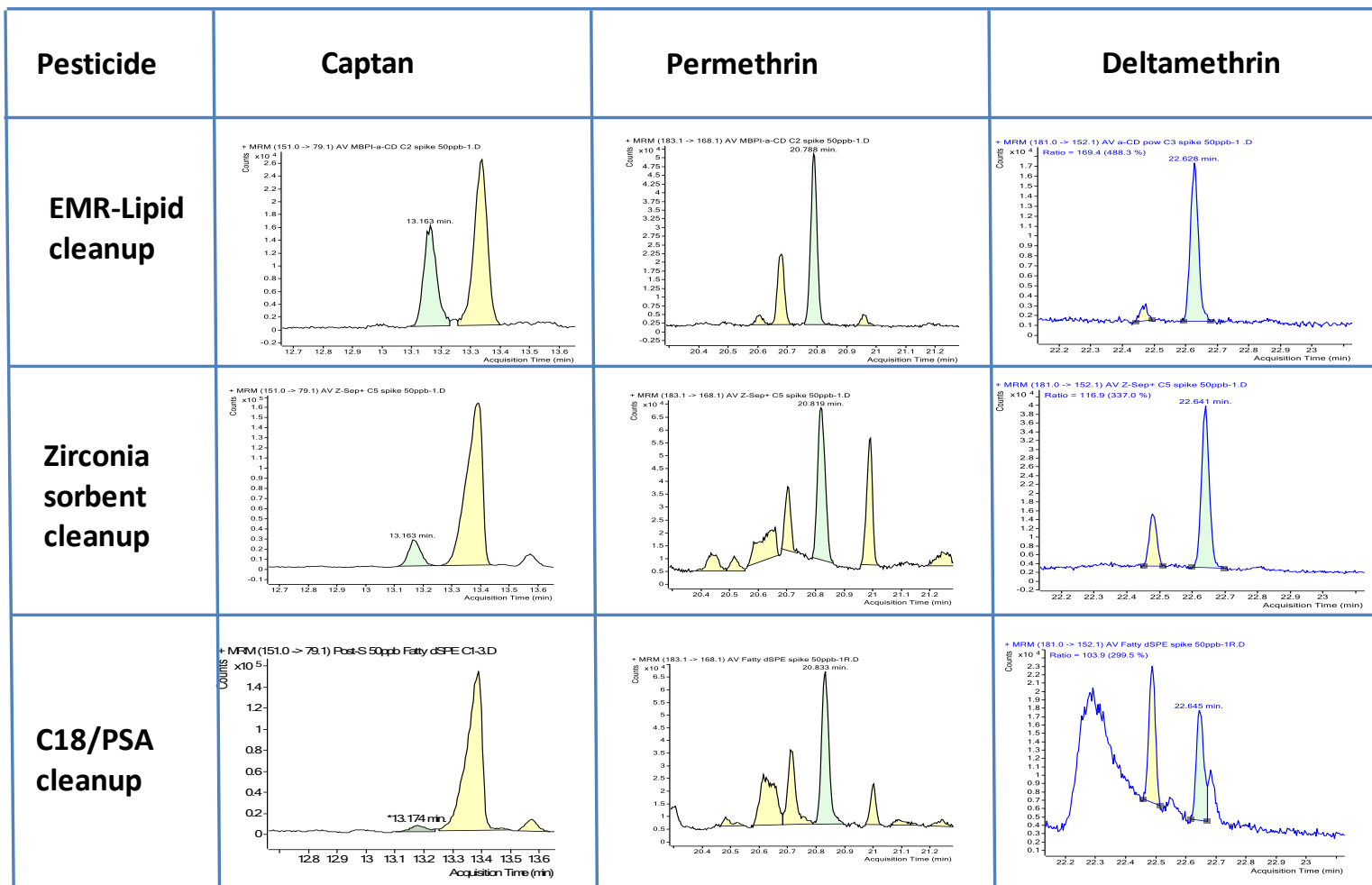

5994-0405EN

# Analysis of Multiclass Multiresidue Pesticides in Milk Using Agilent Captiva EMR-Lipid with LC/MS/MS and GC/MS/MS (5994-2038EN)

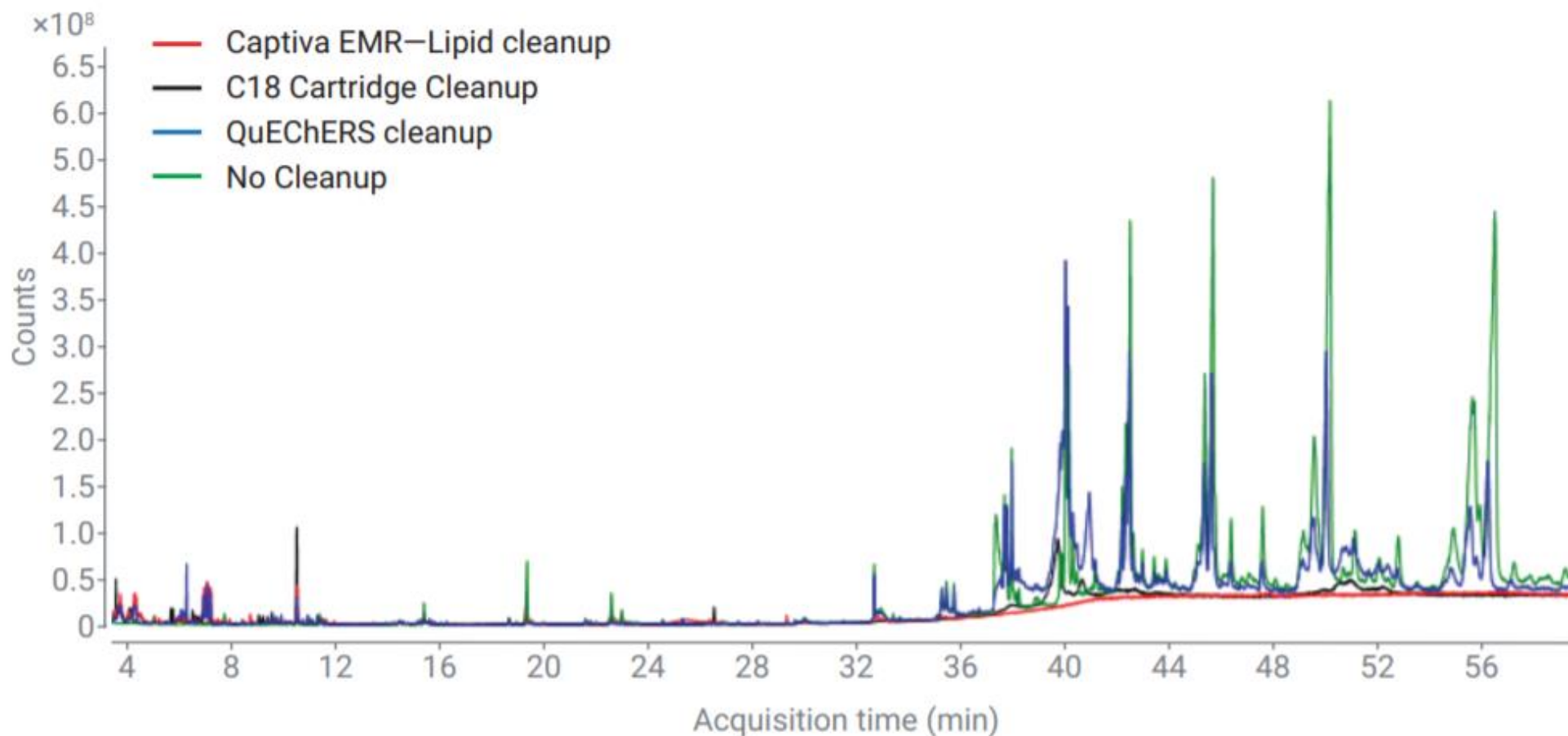

# Carbon Material Used in Food Analysis

- The structure of carbon material, including graphitized carbon black (GCB), coconut carbon, and activated carbon, favors the retention of pigment components.
- Pigments are one of major matrix interferences in plant-origin food matrices
- Offers efficient pigment removal, but can also cause unwanted targets loss, such as planar compounds



Color	Pigments
Green	Chlorophyll Lutein
Red, blue, purple, black	Anthocyanidins Anthocyanins
Orange, yellow	Carotenoids Xanthophylls


# Carbon S Sorbent

Agilent Carbon S sorbent is an *advanced hybrid carbon material* with **optimized carbon content and pore structure**

Significant improvement on sensitive pesticides recovery and reproducibility


Decent pigment removal, comparable to GCB sorbent and better than other sorbents for pigment removal

Improved sorbent selectivity and thus better balance

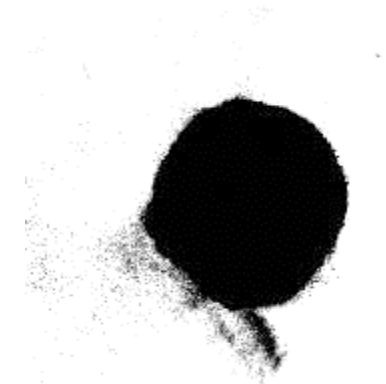


Carbon S sorbent

# Chemical Filtration – Captiva EMR with Carbon S

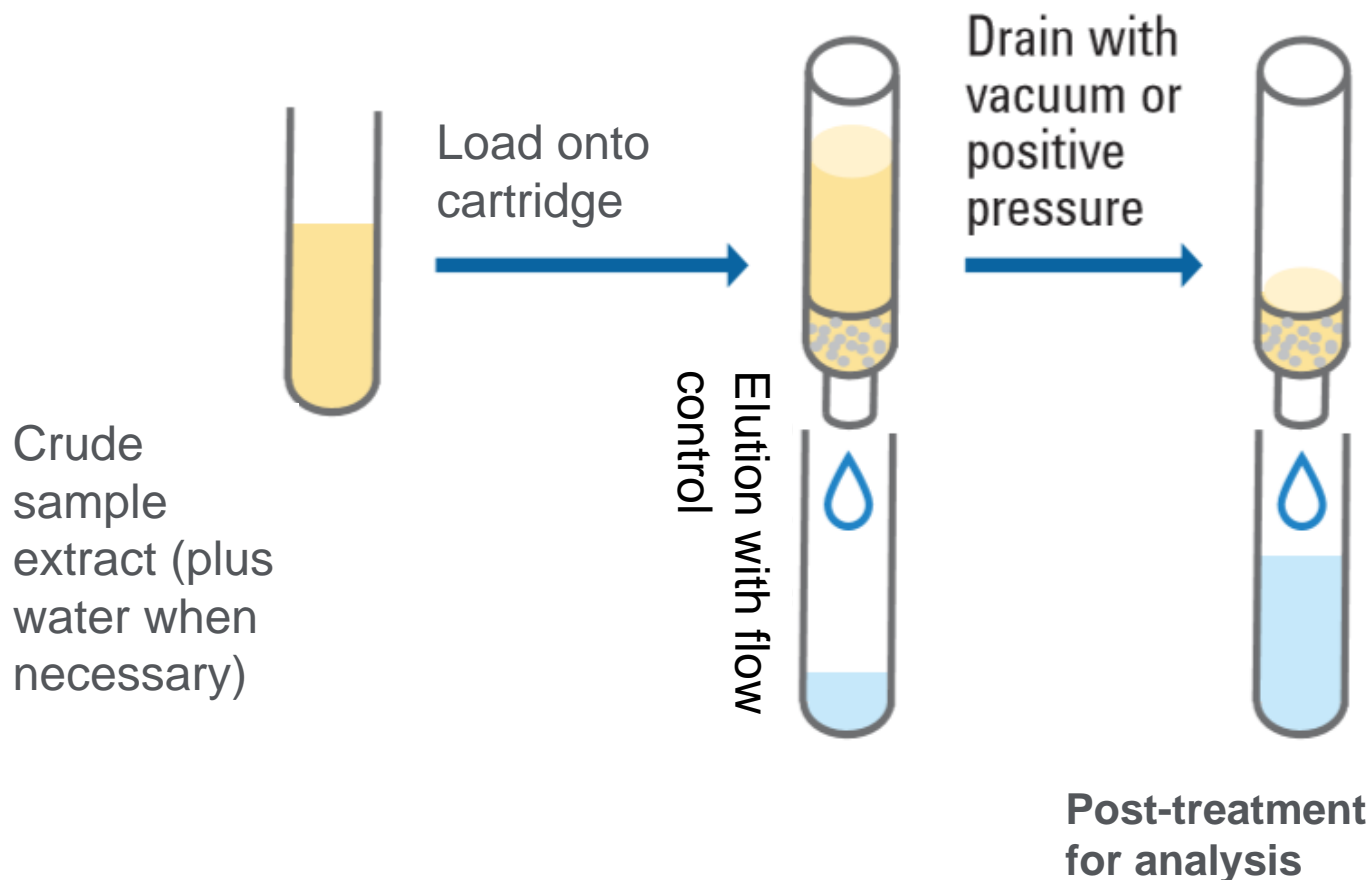

- Sorbent interaction is targeted to unwanted matrix interferences
- Limited interactions with targets to prevent analytes loss
- Basic methodology
  - Load the sample crude extract to cartridge
  - Elution with flow control, either on gravity or low-level external force
  - Collect eluent for direct analysis or further post-treatment

# Captiva EMR with Carbon S Passthrough Cleanup vs dSPE Cleanup

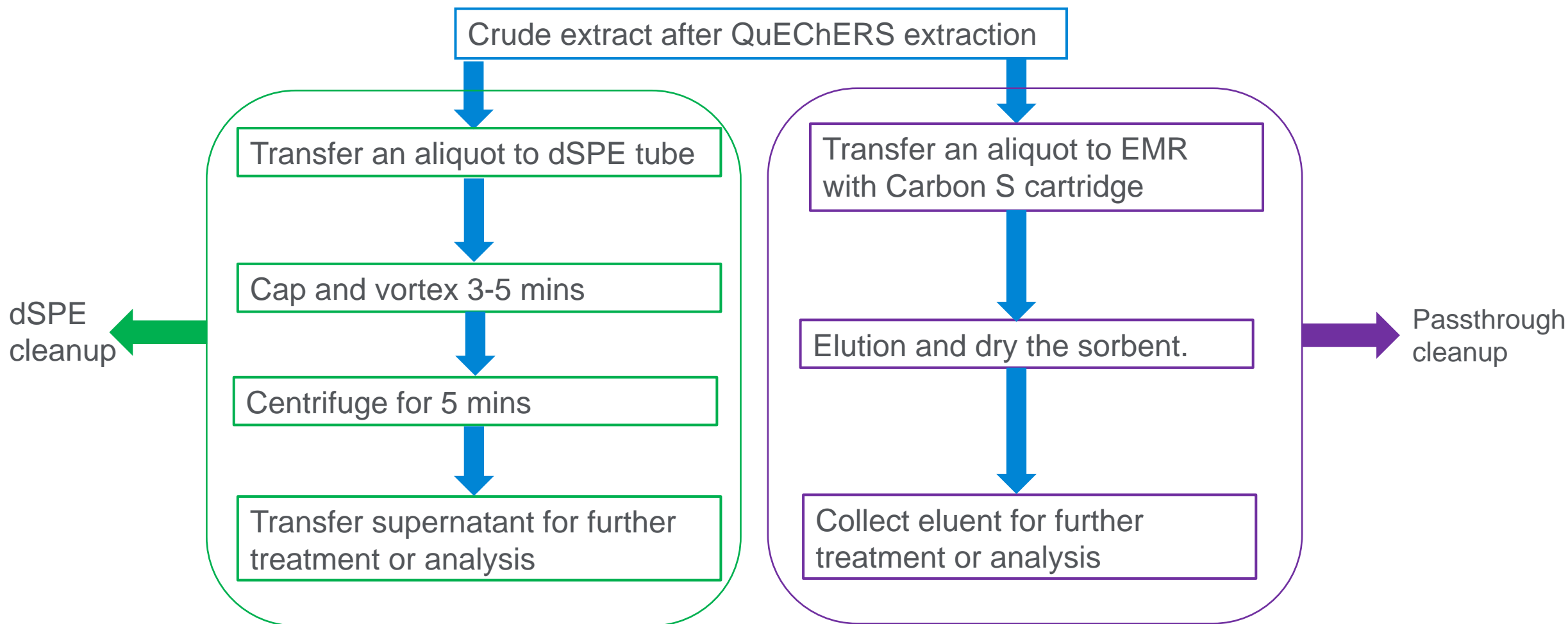

Passthrough cleanup saves 15 to 30% of processing time. The more samples to process, the more time can be saved.

# Captiva EMR Passthrough Selection Guide



**Captiva EMR-Lipid**  
High Lipid/Fats/Oils  
• Meat, dairy, oils, eggs



**Captiva EMR-HCF**  
High Chlorophyll Fresh  
• Spinach, arugula, chard



**Captiva EMR-GPF**  
General Pigmented Fresh  
• Berries, peppers, broccoli




**Captiva EMR-GPD**  
General Pigmented Dry  
• Spices, seasoning, herbal medicine



**Captiva EMR-LPD**  
Low Pigmented Dry  
• Nuts, tobacco, light pigmented spices



# Captiva EMR – Cleaner Matrix Background on GC/MS/MS

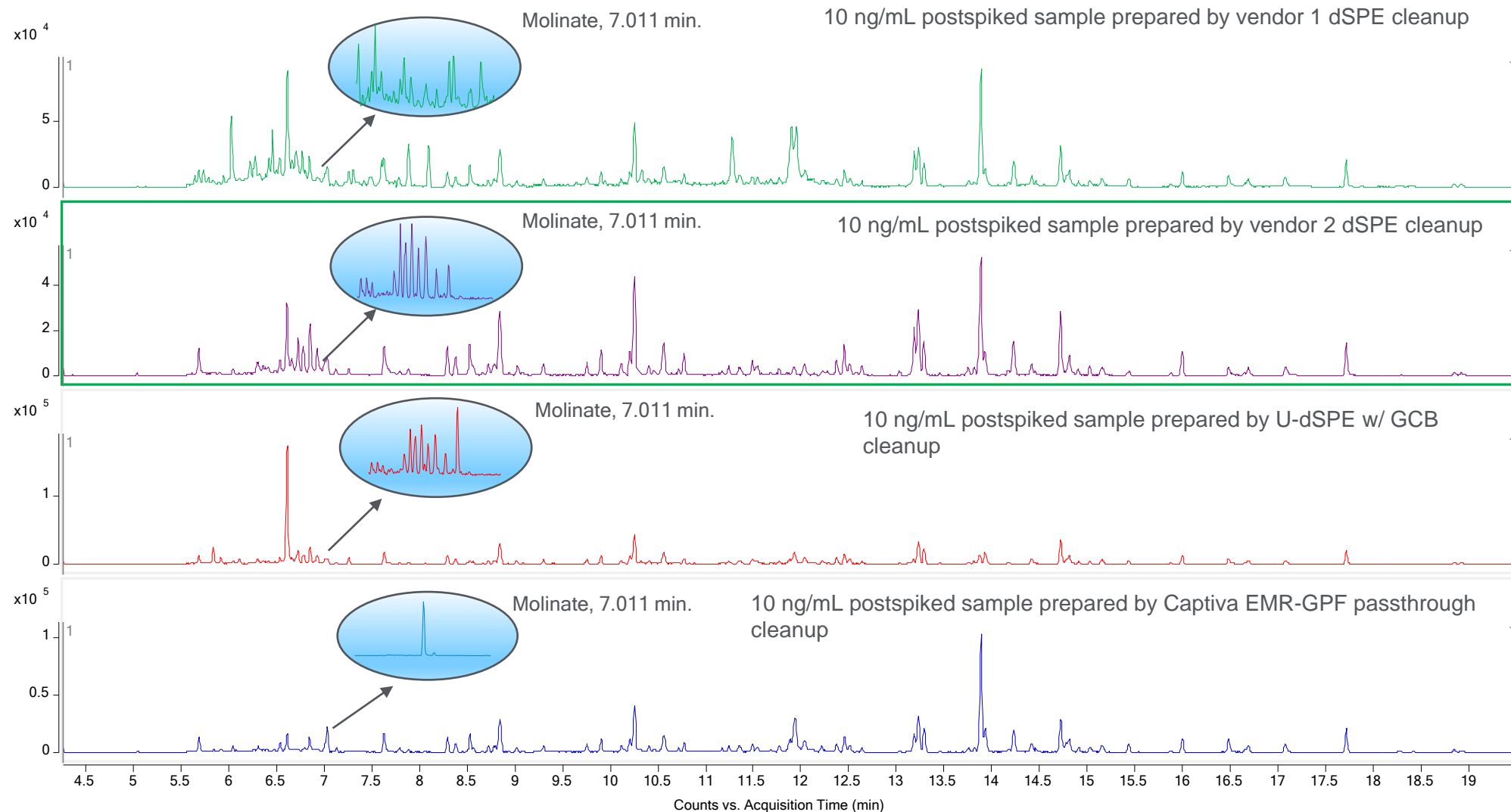



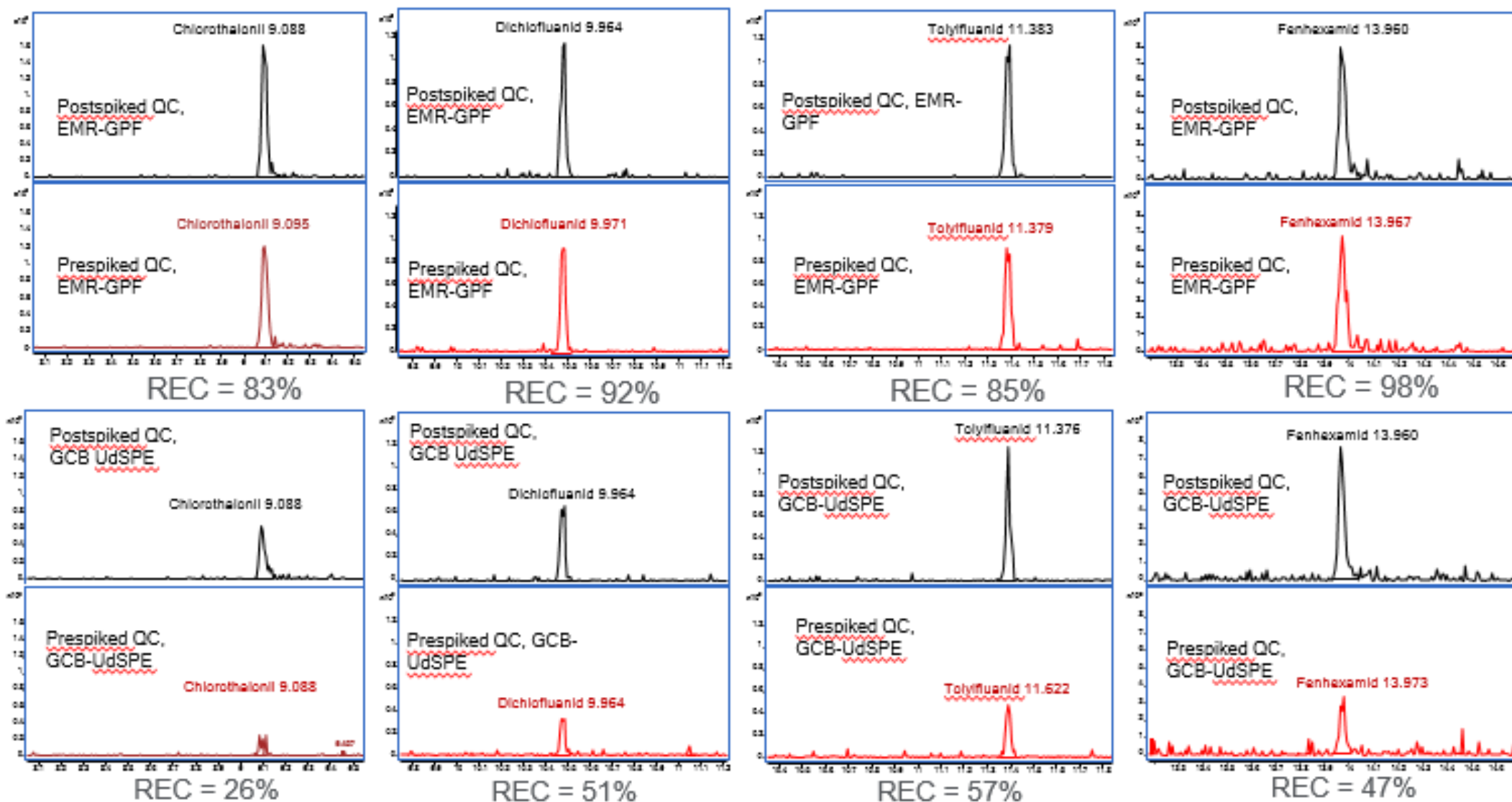


Figure 6. Chromatographic comparison of sensitive targets for samples prepared using different cleanup method.


# Drawbacks of Liquid-Liquid Extraction

- LLE does have drawbacks
  - Inconsistent results from one analyst to another
    - Shaking time
    - Shaking motion
    - Determination of where to cut between layers
  - Emulsions
  - Labor intensive
  - Quite tedious with small sample sizes (<5 mL)
  - Challenging with large numbers of samples
  - Difficult to automate for large numbers of samples

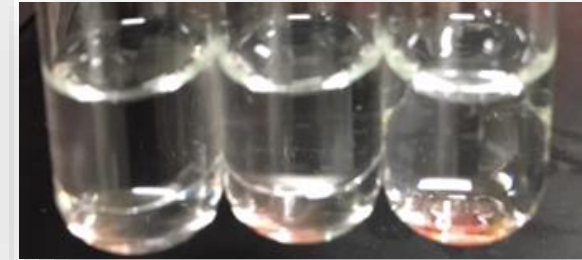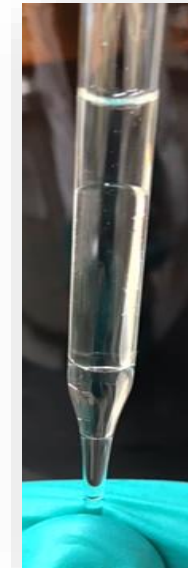

How many of these problems can be fixed with solid supported liquid extraction?


# How Does SLE Work?

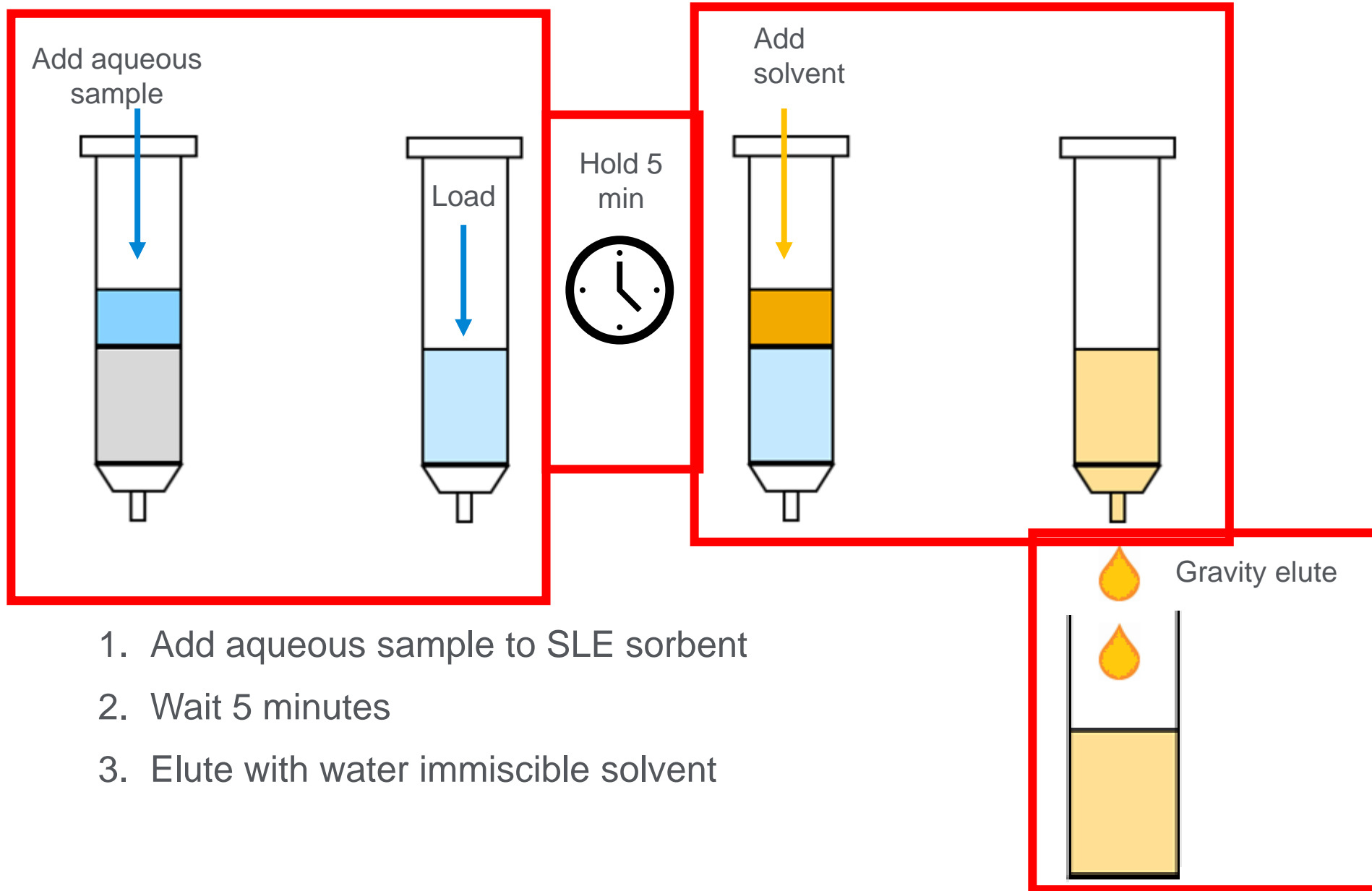

1. Add aqueous sample to SLE sorbent
2. Wait 5 minutes
3. Elute with water immiscible solvent

# What is SLE Sorbent?


- There are two types of SLE media
  - Diatomaceous earth (DE) based products like our Chem Elut brand of SLE products
    - A mined fossil diatom material, which is heterogeneous and inconsistent from one mine to the next

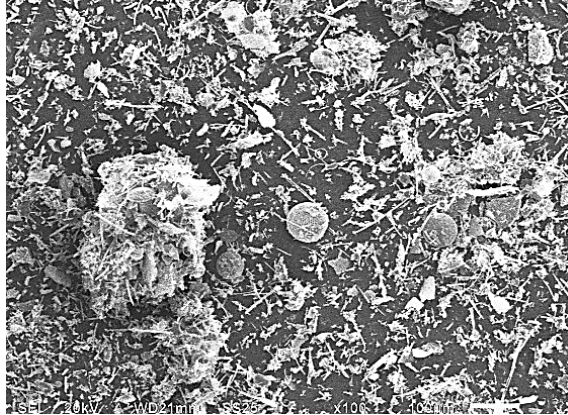

## Diatomaceous earth in Chem Elut

- ✗ Naturally occurring; mined
- ✗ Broad particle size distribution
- ✗ Supplier reliability issues
- ✗ Poor lot-to-lot consistency

- Synthetic media we use in Chem Elut S
  - Controlled synthesis to be consistent batch after batch

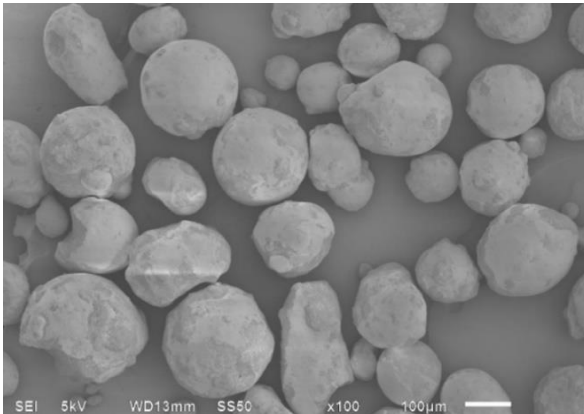

## Synthetic SLE sorbent

- ✓ Large scale synthesis
- ✓ Narrow particle size distribution
- ✓ Reliable supplier
- ✓ Controlled manufacturing

# Supported Liquid Extraction (SLE)

## Chem Elut S

- Same extraction mechanism as traditional liquid-liquid extraction (LLE)
- Cartridge and plate format, packed with proprietary synthetic sorbent – high surface area
- Simple method, gravity flow
- Smaller volume sample and solvent compared to LLE
- No emulsions

Cartridges for sample volumes 0.2 – 20 mL



Bulk Chem Elut S  
1 and 4 kg


96-well plate for sample  
volume 200 and 400  $\mu$ L

# Solid Phase Extraction (SPE)

- Capabilities
  - Very selective
  - Highly clean samples
  - Concentrated samples
  - Wide range of applicability
  - Automation friendly
- Types of SPE
  - Nonpolar (reversed phase) SPE
  - Polar (normal phase) SPE
  - Cation exchange SPE
  - Anion exchange SPE
  - Mixed mode SPE
  - Specialty SPE

Bond Elut:

Silica or polymer based, cartridge and 96-well plate format

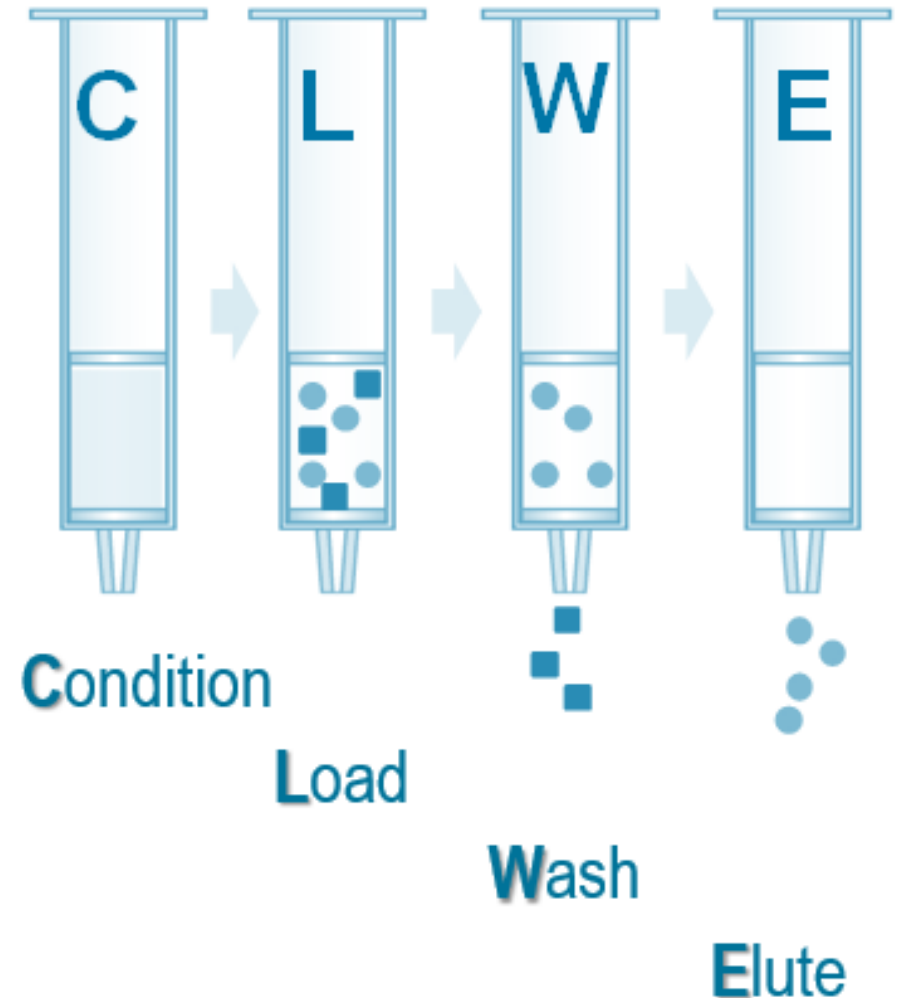

# Agilent SPE Offering

- Reliable SPE with a 30-year history
- Agilent offers the most comprehensive set of phases, sizes, and formats of any SPE provider (over 40 sorbent materials/phases available).
- Easy adoption of methods due to high number of publications and applications
- Includes packed bed silica and polymeric phases, and monolithic silica phases

## **OMIX** monolithic silica tip SPE

OMIX C18  
OMIX MP1  
OMIX SCX

## **SPEC** monolithic silica disk SPE

SPEC C2  
SPEC C8  
SPEC C18  
SPEC C18AR  
SPEC PH  
SPEC NH2  
SPEC CN  
SPEC Si  
SPEC PSA  
SPEC SAX  
SPEC SCX  
SPEC MP1  
SPEC MP3

## **Bond Elut Silica and polymer SPE**

Bond Elut AccuCAT  
Bond Elut Alumina (AL-A)  
Bond Elut Alumina (AL-B)  
Bond Elut Alumina (AL-N)  
Bond Elut NH<sub>2</sub>  
Bond Elut C1  
Bond Elut C2  
Bond Elut C8  
Bond Elut C18 .....  
..... **40 phases**

## **Bond Elut Plexa polymer SPE**

Bond Elut Plexa  
Bond Elut Plexa PCX  
Bond Elut Plexa PAX




# SPME Fiber and Arrow Offering from Agilent

## Solid phase microextraction (SPME)

- Environmental analyses of water samples
- Odor analyses (ppt)
- Flavor analyses of food products
- Surfactants, other industrial applications
- Trace analysis in food
- Herbicides/pesticides
- Trace impurities in polymers and solid samples
- Solvent residues in raw materials
- Explosives



SPME fibers

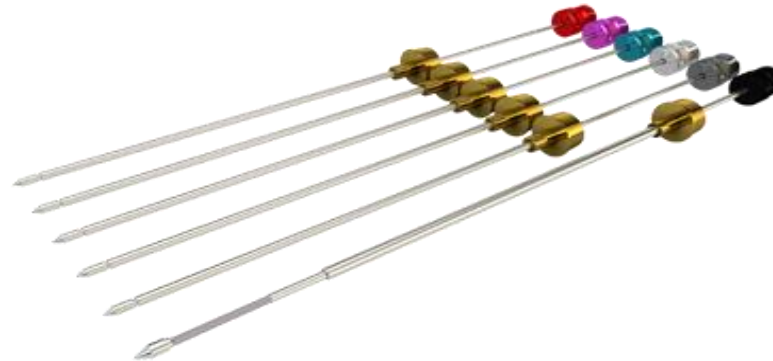


SPME arrows

Smart SPME Fibers and Arrows are now available for the smart rail systems


# Examination of Lower Molecular Weight PAHs in Drinking Water Using Agilent PDMS SPME Fibers

Polycyclic aromatic hydrocarbons (PAHs) are a large class of organic compounds containing two or more fused aromatic rings. PAHs are considered compounds of concern by environmental organizations; their concentration in water is strictly regulated.

5994-1301EN

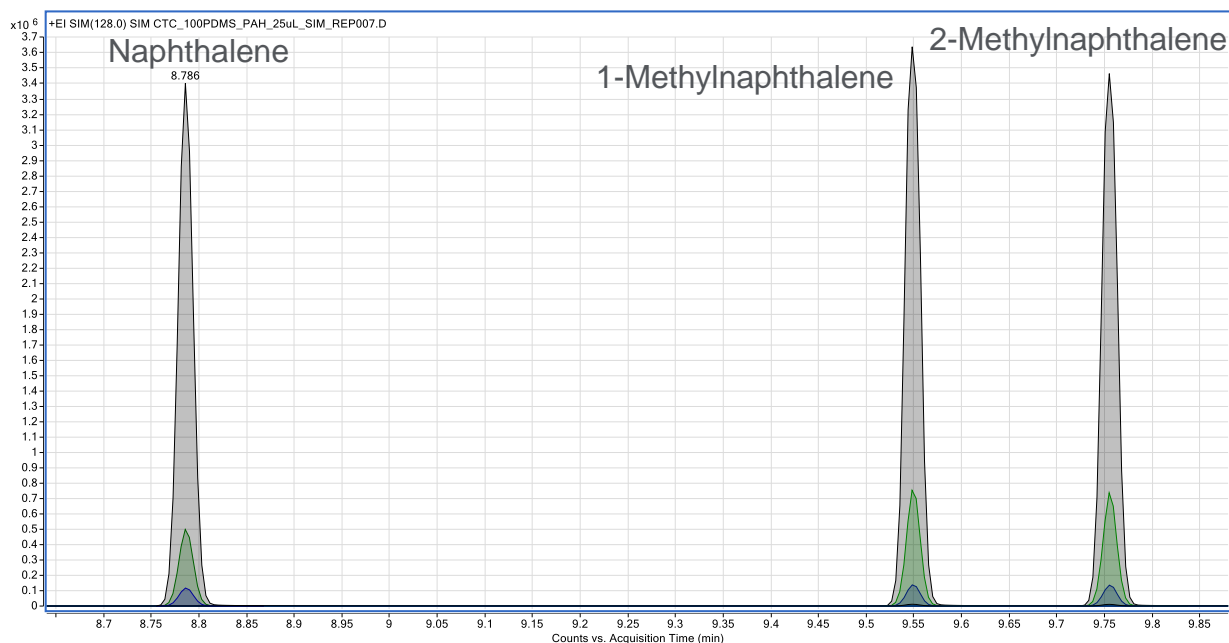

SIM chromatogram of naphthalenes with PDMS fibers (black trace = 100 µm; green trace = 30 µm; blue trace = 7 µm)

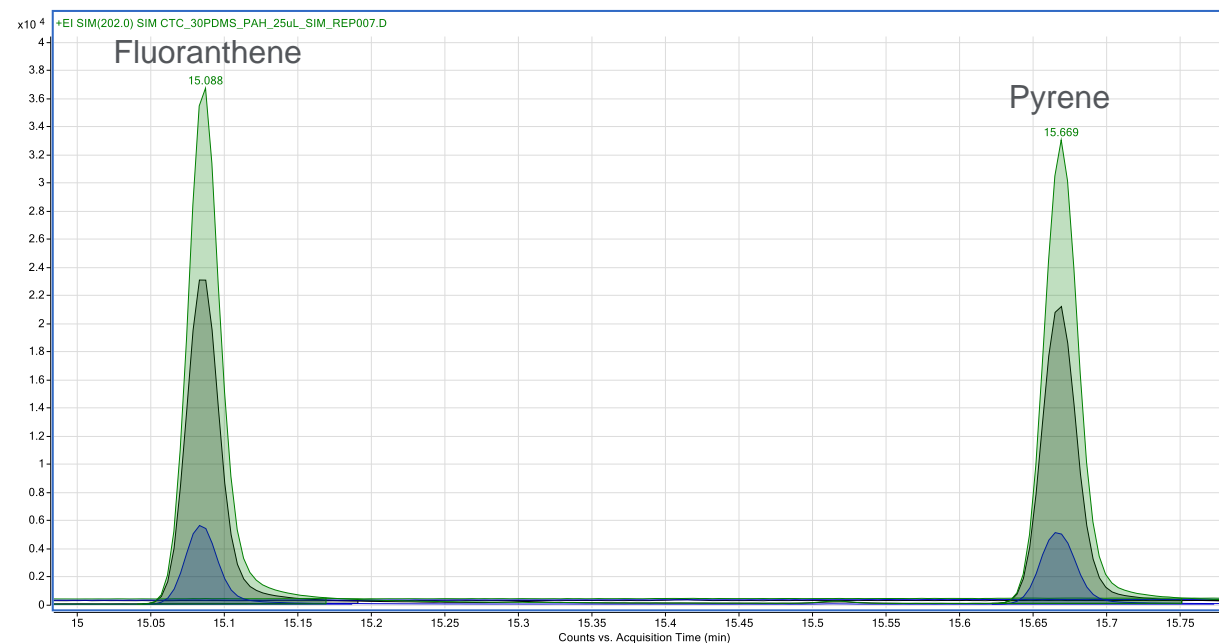

SIM chromatogram of fluoranthene and pyrene with PDMS fibers (black trace = 100 µm; green trace = 30 µm; blue trace = 7 µm)

# Agilent Bond Elut QuEChERS

Quick Easy Cheap Effective Rugged and Safe

Initially developed for screening of pesticide residues in fruit and vegetables to make sample cleanup of food faster, simpler, less expensive, and greener.

Now, QuEChERS is used with other matrices and compound classes as well.

Consists of two steps, and therefore two kits:

Step 1: Liquid extraction



Step 2: Dispersive SPE/  
interference removal



# QuEChERS Workflow

QuEChERS extraction salts

QuEChERS dispersive  
SPE sorbents

## Step 1: Salting Out Extraction



1 Weigh sample



2 Add water and DC spikes if needed and spike with internal standard



3 Add acetonitrile



4 Vortex or shake



5 Add salt packet



6 Shake 1 minute



7 Centrifuge at 4000 rpm for 5 minutes



Phase separation of acetonitrile and aqueous layer

## Step 2: Dispersive Solid Phase Extraction (dSPE)



1 Choose the dispersive cleanup kit and add acetonitrile extract



2 Vortex for 1 minute



3 Centrifuge at 4000 rpm for 5 minutes



4 Take aliquot of supernatant and dry down or dilute as necessary



5 Place in autosampler vials for GC or LC analysis

## Step 3: Analysis

Analysis GC or LC MS, MS/MS

# Bond Elut Dispersive SPE Kits



## Dispersive kit

Centrifuge tubes containing preweighed SPE sorbent such as:

- C18: Removes residual fats and lipids
- PSA: 'Primary/secondary amine' for removal of organic acids and sugars
- Carbon S: Graphitized carbon black, removes pigments
- EMR-Lipid: Removes unbranched hydrocarbon chains (lipids)

Dispersive SPE kits are available for different food types.

They are for both AOAC (US) method and EN (Europe).

QuEChERS is a nonselective technique and does not remove **all** matrix, just enough.

Dispersive sorbents are also available as bulk material.

# Sample Preparation Application Note Examples using GC Detection

## Captiva EMR-Lipid

1. [Determination of 14 Polycyclic Aromatic Hydrocarbon Compounds in Edible Oil, PN 5994-1483EN](#)
2. [Analysis of Multiclass Multiresidue Pesticides in Milk, PN 5994-2038EN](#)

## Captiva EMR with Carbon S

1. [Determination of Over 300 Pesticides in Cayenne Pepper, PN 5994-5630EN](#)
2. [Determination of Multiclass, Multiresidue Pesticides in Bell Peppers, PN 5994-4767EN](#)

## Chem Elut S

1. [Determination of Aromatic Amines Derived from Azo Colorants by GC/MS, PN 5994-0951EN](#)

## SPME

1. [Analysis of Free Volatile Phenols in Smoke-Impacted Wines by SPME, PN 5994-3161EN](#)



# Manifolds for Processing Cartridges and 96-Well Plates

Captiva vacuum collar



Positive Pressure Manifolds

SPS 24 vacuum manifold



Vac Elut 20 vacuum manifold




Vac Elut 12 vacuum manifold



96 well plate vacuum manifold



# Guard Column or Retention Gap

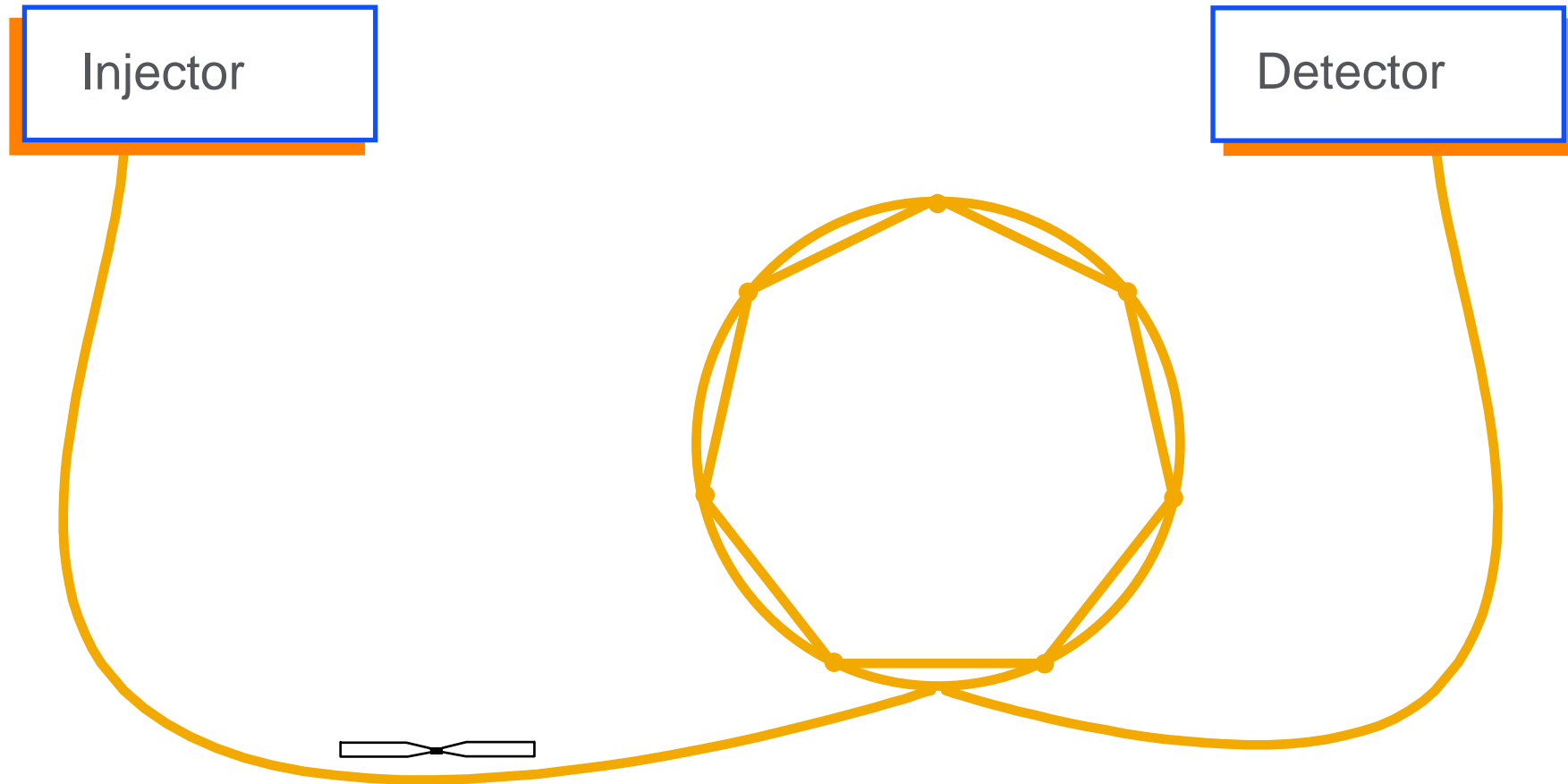


The guard column is 3 to 5 m of deactivated fused silica tubing with the same diameter as the analytical column. It is connected with a zero dead volume union.



# Better Connections: Ultra Inert Press Fits or Ultimate Union

- Ultra Inert press fits:
  - Join retention gap or guard column to analytical column
  - Dependable inertness performance at a lower cost
  - Batch certified inertness
  - Improved packaging and installation instructions
  - Easier to use – transparent deactivation gives visibility of the column connection
- Ultimate union
  - More robust
  - Reusable
  - Recommended for users with MS

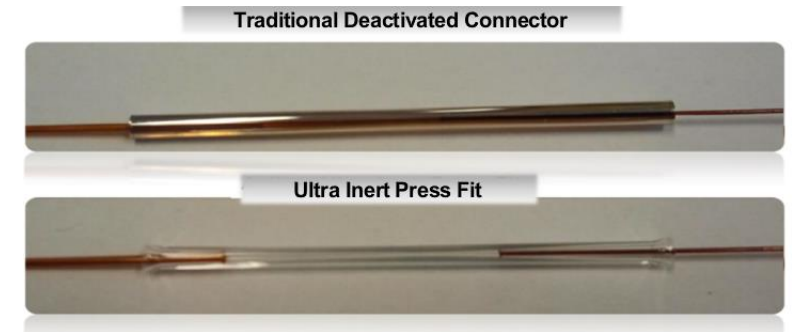

# Integrated Guards – DuraGuard

- No union
- Possible for any DB column 0.18 mm and larger
- Limited offering “off-the-shelf”

## DuraGuard


Phase	ID (mm)	Length (m)	Film (µm)	Guard Length (m)	Part No.
DB-1	0.25	30	0.25	10	122-1032G
DB-XLB	0.25	30	0.25	10	122-1232G
DB-5ms	0.25	30	0.25	10	122-5532G
			0.50	10	122-5536G
			1.00	10	122-5533G
		60	0.25	10	122-5562G
	0.32	30	1.00	10	123-5533G
	0.53	30	0.50	10	125-5537G
DB-5.625	0.18	20	0.36	5	121-5622G5
	0.25	30	0.25	5	122-5631G5
DB-1701	0.53	30	1.00	10	125-0732G
DB-624	0.53	30	3.00	5	125-1334G5

# An Introduction to Backflush Techniques




# Why Use Backflush?

- Avoid unwanted sample components from entering analytical column
- Avoid heavy compounds from reaching detector
- Shorten run times and increase sample throughput

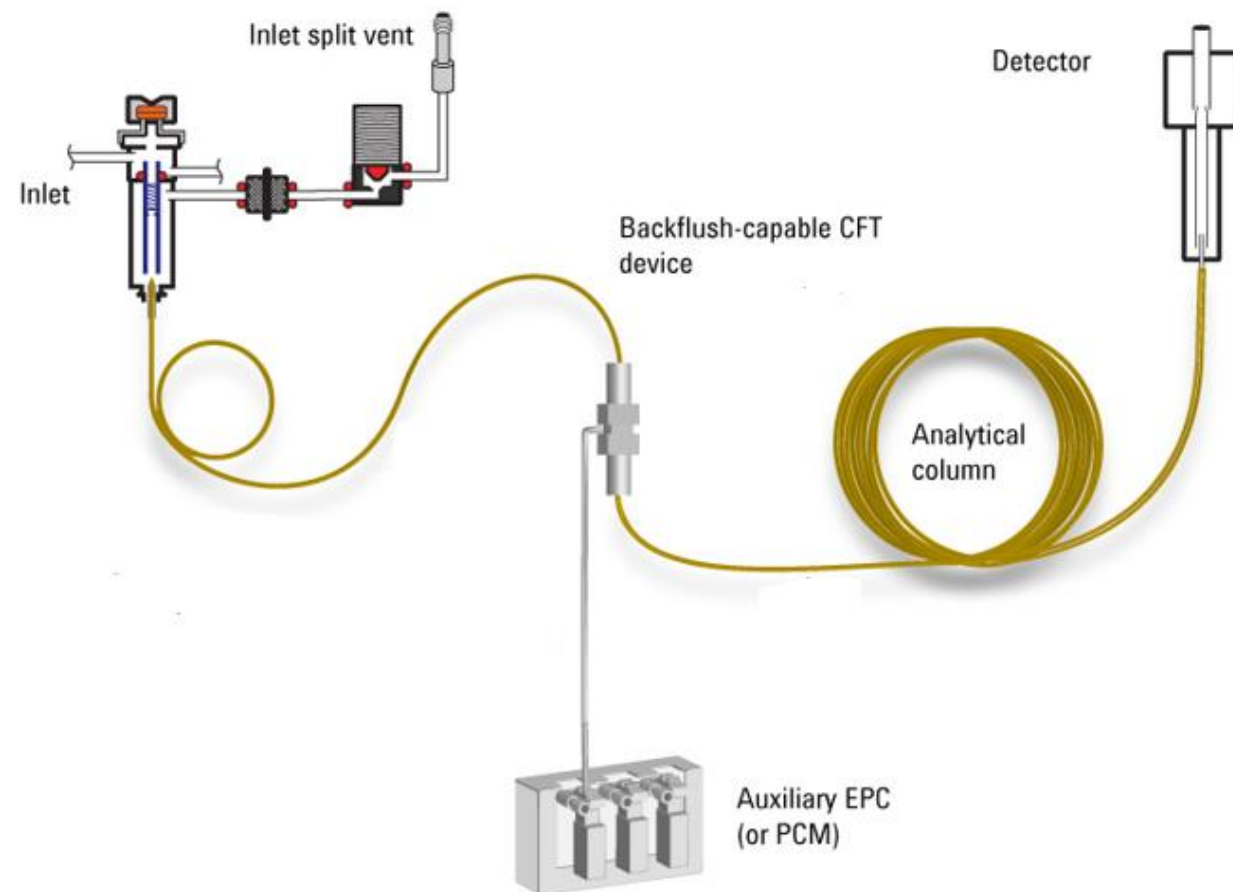

# Why Use Backflush?


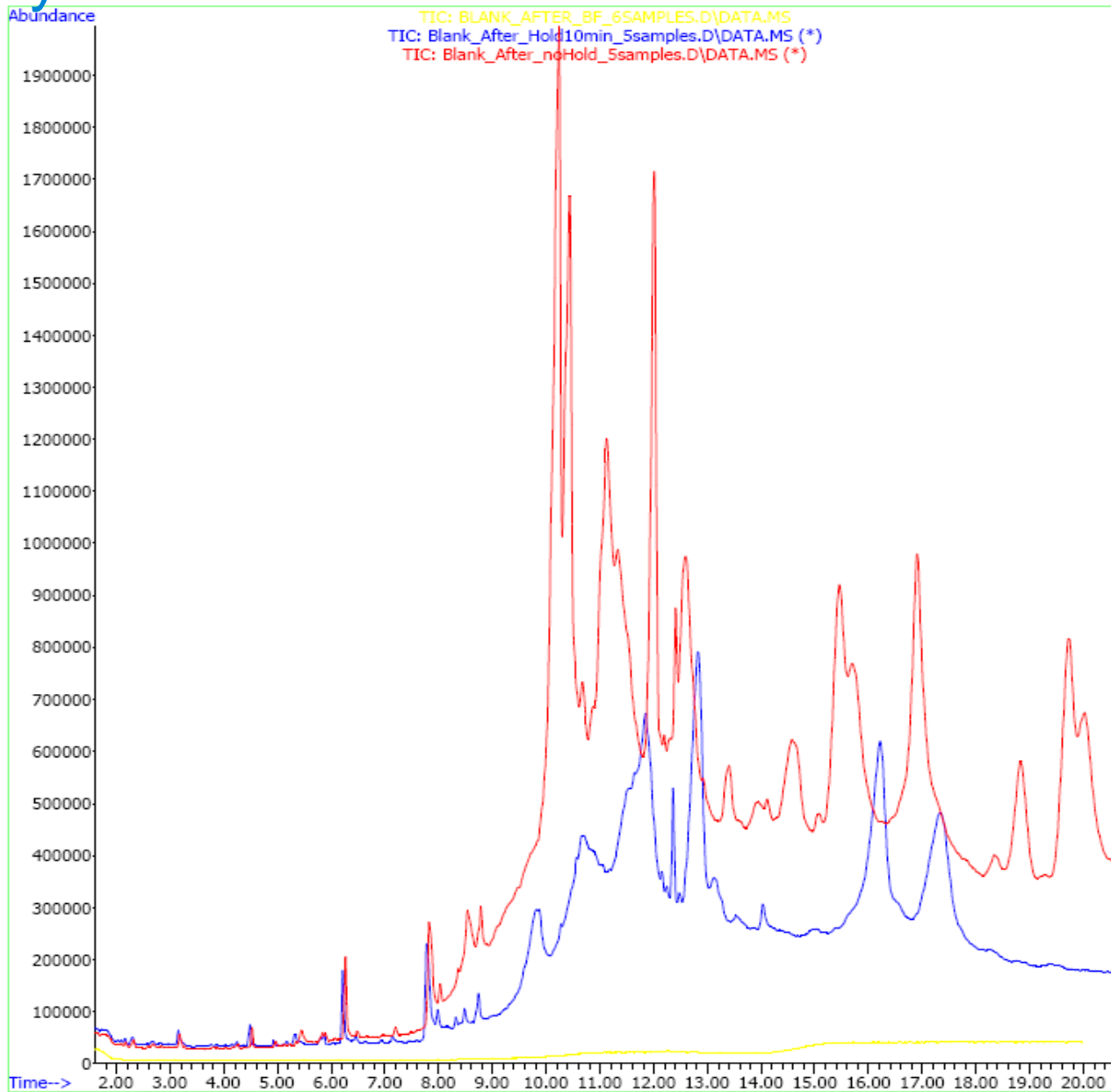


**Blank after five injections  
(normal method)**


**Blank after five injections  
(with 10 minute bake at 280 °C)**

**Blank after five injections  
(2 minute backflush, 50 PSI at 280 °C)**


# Why Use Backflush?

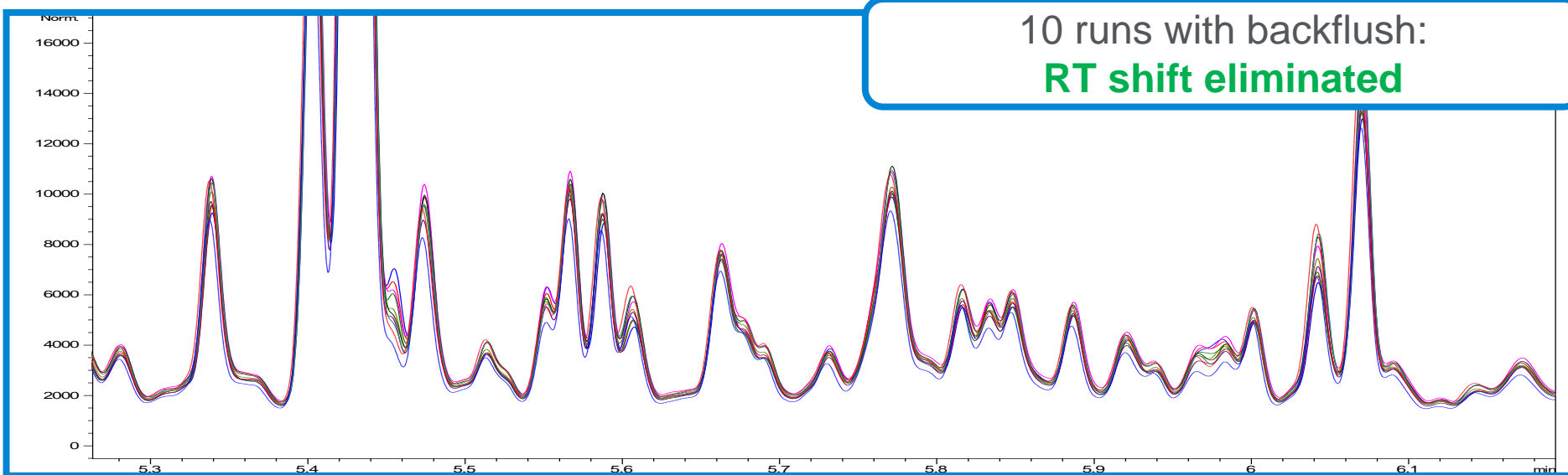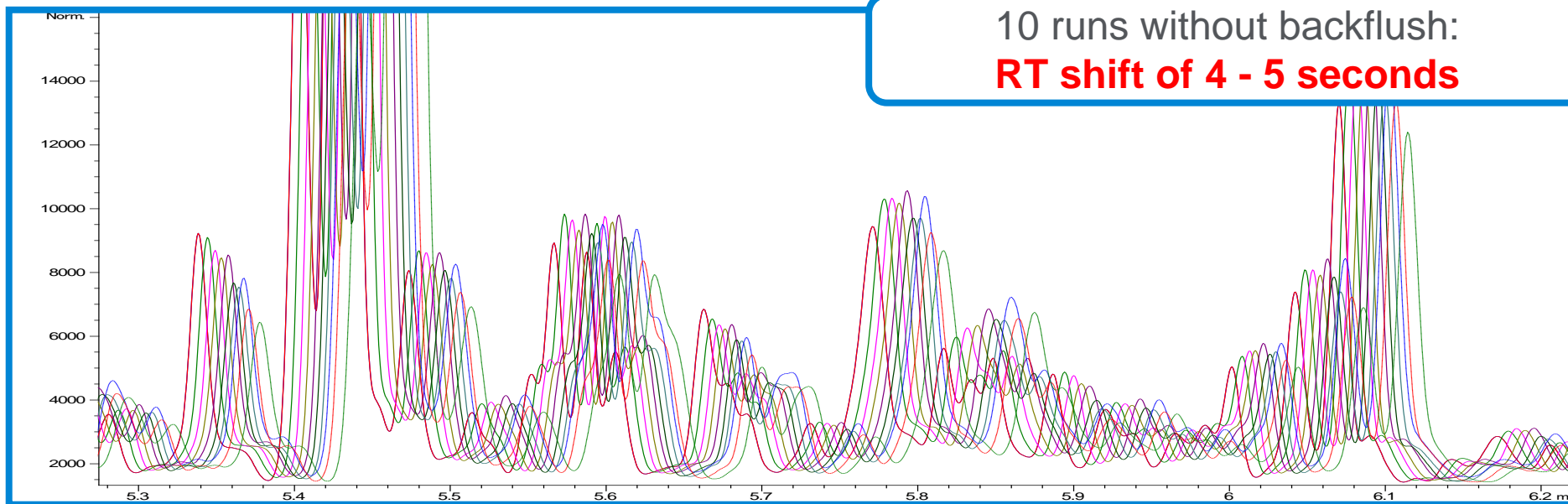

# Why Use Backflush?

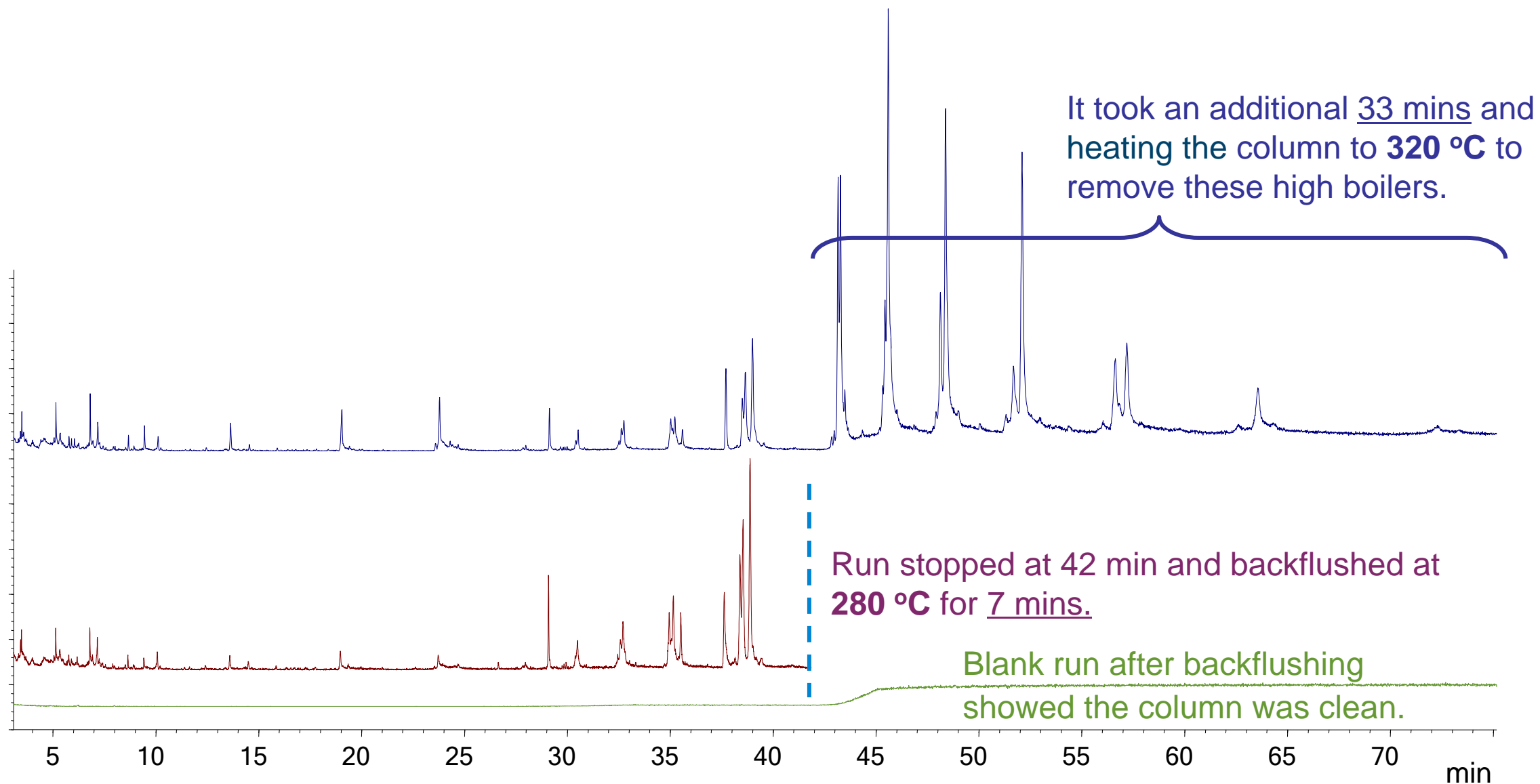

# Why Use Backflush?

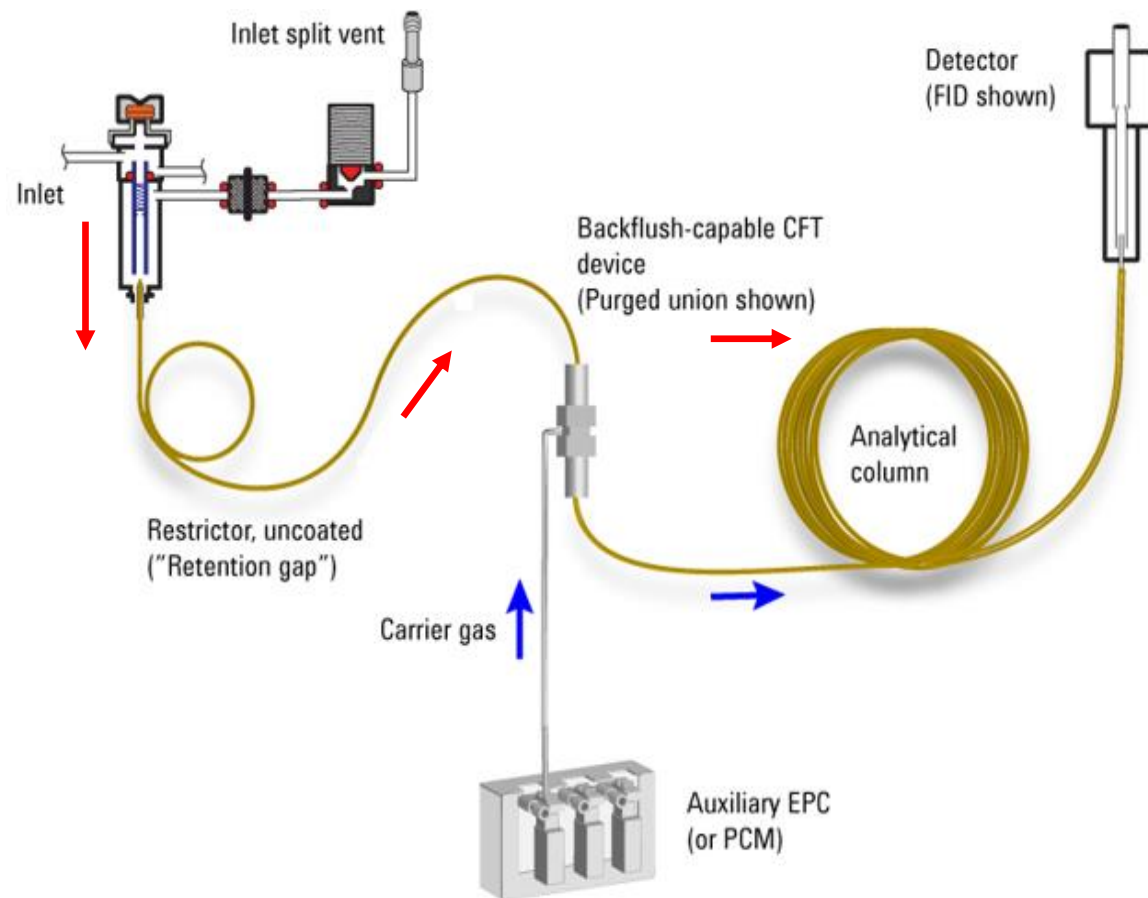



# Summary

- Stay within the temperature limits of your column and handle it carefully
- Be mindful of what you are injecting into your system
- Take notice of any chromatography or baseline changes from when the column was brand new
- Ultimately, sample preparation/cleanup is the most reliable way to address common chromatography data problems.
- Agilent offers a wide range of sample preparation products to support your analysis using established methods and protocols:
  - Filtration, protein and lipid removal
  - SLE
  - QuEChERS
  - SPE
- Backflush techniques can help extend the lifetime of your column if your instrument can be set up with it

# Contact Agilent Chemistries and Supplies Technical Support



1-800-227-9770 option 3, option 3:

Option 1 for GC and GC/MS columns and supplies

Option 2 for LC and LC/MS columns and supplies

Option 3 for sample preparation, filtration, and QuEChERS

Option 4 for spectroscopy supplies

Option 5 for chemical standards

**Available in the U.S. and Canada 8–5, all time zones.**



[gc-column-support@agilent.com](mailto:gc-column-support@agilent.com)

[lc-column-support@agilent.com](mailto:lc-column-support@agilent.com)

[spp-support@agilent.com](mailto:spp-support@agilent.com)

[spectro-supplies-support@agilent.com](mailto:spectro-supplies-support@agilent.com)

[chem-standards-support@agilent.com](mailto:chem-standards-support@agilent.com)