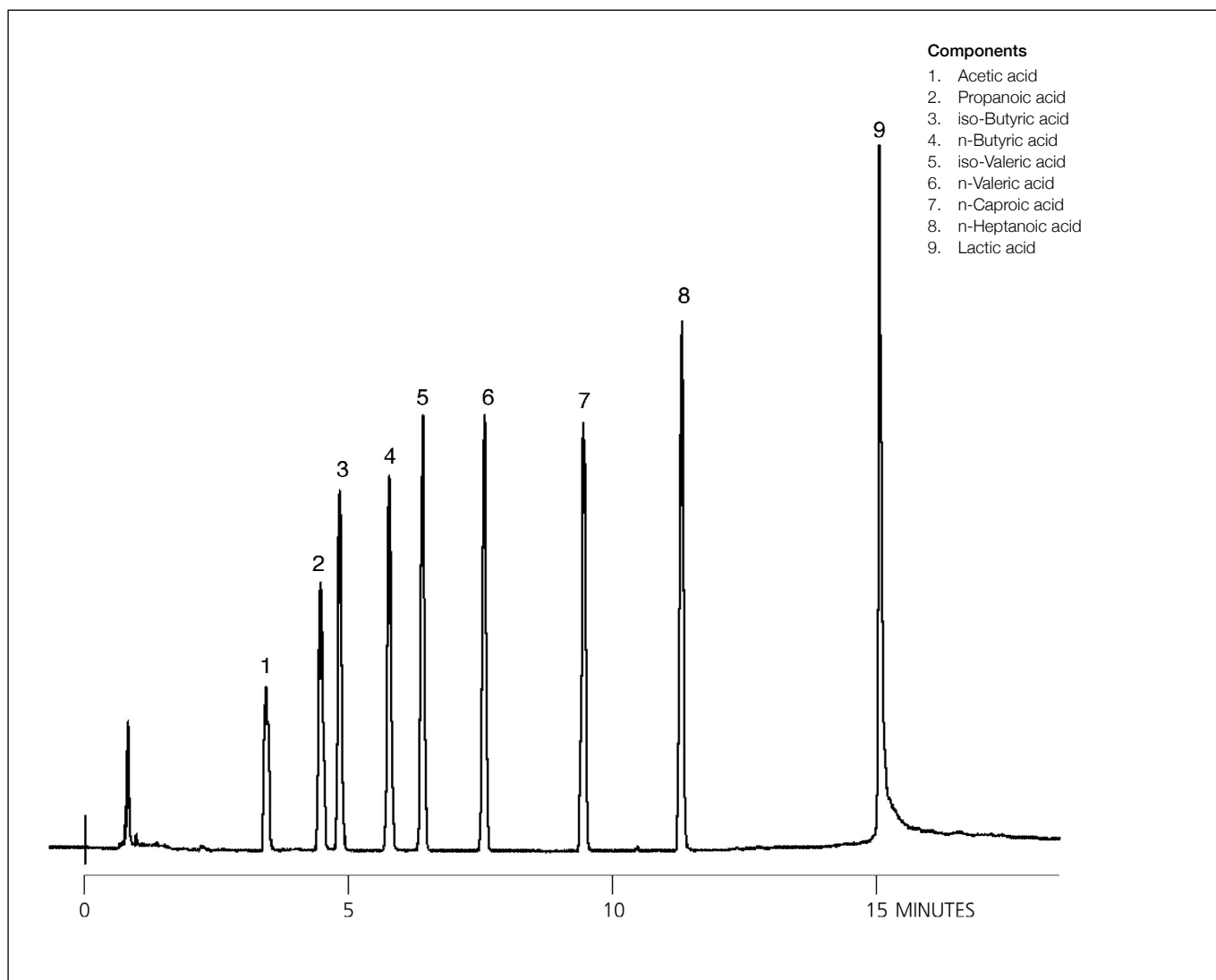


Analysis of organic acids in water

BP21 (FFAP)

Column part number	054477		
Phase	BP21 (FFAP)	Rate	6°C/min
Column	30 m x 0.53 mm x 0.5 µm	Final temperature	180°C, 5 min
Sample	Organic acids in water (with 0.03 M Oxalic acid)	Detector	FID
		Injection mode	On-column
Initial temperature	85°C, 0 min	Carrier gas	H ₂


 For more information visit www.trajanscimed.com or contact techsupport@trajanscimed.com