



Drugs of abuse

Application Note

Forensic Toxicology

Authors

Agilent Technologies, Inc.

Introduction

Fast GC separation of four drugs of abuse on an Agilent VF-5ms column in less than four minutes.





Agilent Technologies

Conditions

Technique : GC
Column : Agilent VF-5ms, 0.1 mm x 10 m fused silica
(df = 0.4 μ m) (Part no. CP8934)
Temperature : 140 °C, 0.5 min \rightarrow 300 °C, 30 °C/min
Carrier Gas : Helium, constant flow. Approx. 400 kPa, 0.7 ml/min
Injector : Splitless, 30 s, Temp: 280 °C
Detector : TOF
Interface Temp : 320 °C
Sample Size : 0.5 μ L
Concentration Range : 1 ng/ μ L

Peak identification

1. amphetamine
2. methamphetamine
3. methadone

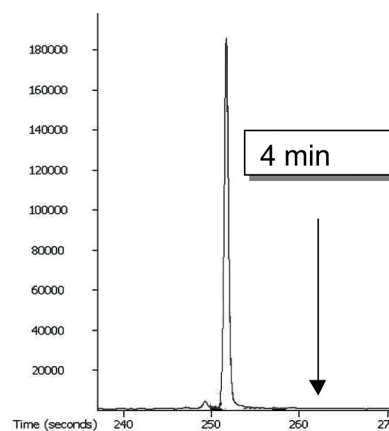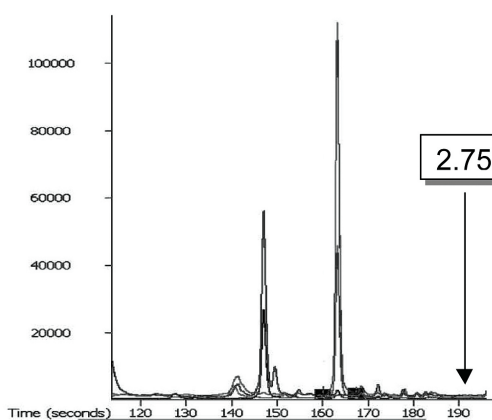

www.agilent.com/chem

For Forensic Use only.

This information is subject to change without notice.

© Agilent Technologies, Inc. 2011

Printed in the USA

31 October, 2011

First published prior to 11 May, 2010

A02070



Agilent Technologies