


InertSearch™ for GC

InertCap® Applications

Pesticides in golf courses

Data No. GA008-0000

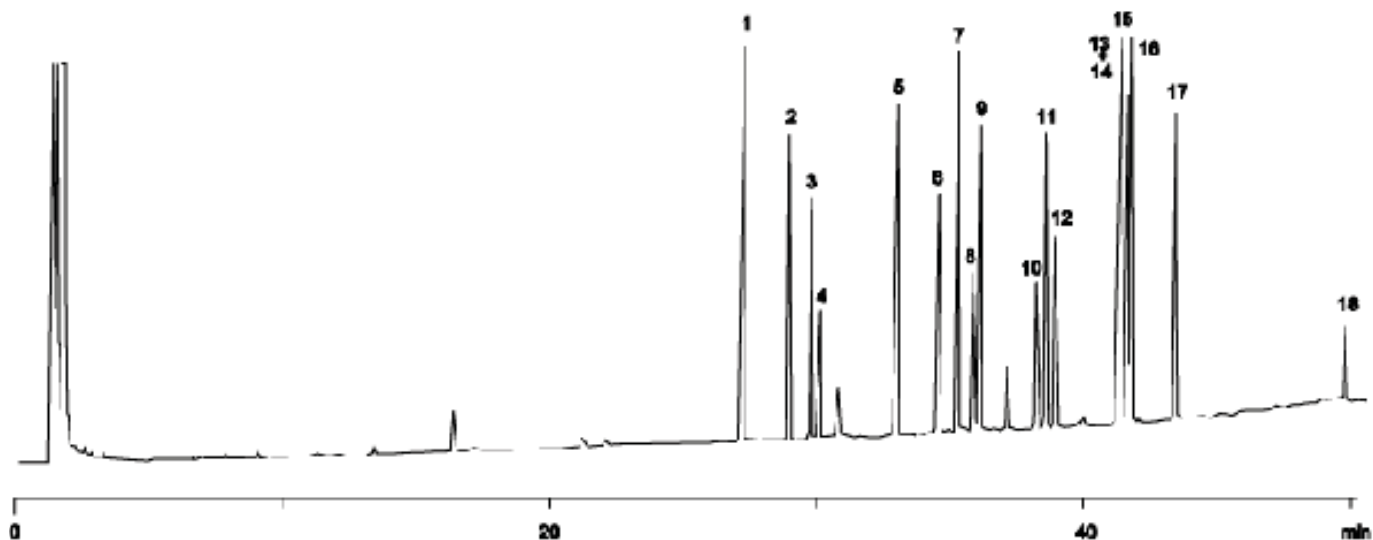

System : GC/FID
Column : InertCap 5
0.25mm I.D. x 30m df=0.4µm
Column Temp. : 100° C → 3° C/min → 250° C
Carrier Gas : He 30cm/sec
Injection : Split 1 : 70
Detection : FID Range 1x8
Sample : Pesticides

Solutes : 1. Simazine(CAT)
2. Propyzamide
3. Diazinon
4. Chlorothalonil(TPN)
5. Tolclophosmethyl
6. Fenitrothion(MEP)
7. Thiobencarb
8. Fenthion(MPP)
9. Chlorpyrifos
10. Pendimethalin
11. Captan
12. Isofenphos
13. Napropamide
14. Butamifos
15. Flutoluanil
16. Isoprothiolane
17. Isoxathion
18. Iprodione