


Improving Data Reliability


LabSolutions DI Assist Package™

Ensuring reliable data and preventing issues such as data tampering are essential requirements in laboratories in a variety of industries. The easy-to-use LabSolutions DI Assist Package extends the know-how cultivated from supporting data integrity, a requirement of the pharmaceutical industry, to instruments from various companies.

The reliability of test data obtained from instruments from various companies is aided by **three automation functions** on a single PC.

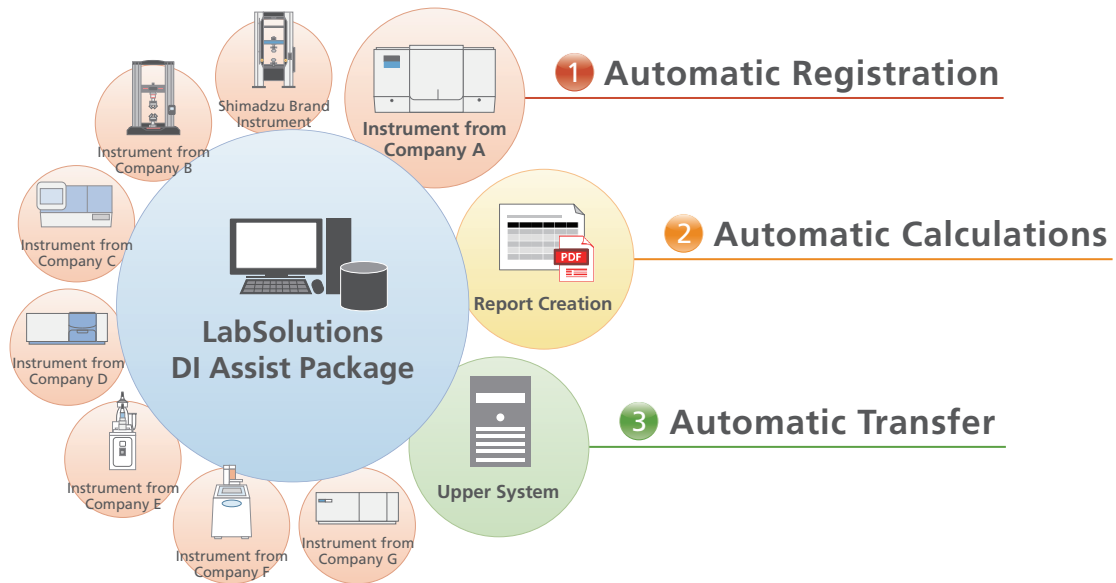


1 Automatic Registration of Test Results Data and PDF Reports from Up to Eight Instruments from Various Companies

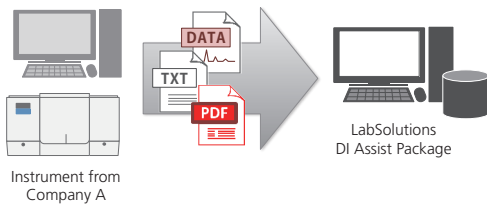

- ▶ Prevents operators from only saving favorable data
- ▶ Makes manual editing during the process impossible
- ▶ Prevents data deletion and concealment


Note 1) This applies only to instruments under PC control.
Note 2) This depends on the specifications for the software from each company, so prior investigation is required.

2 Automatic Calculations by Selecting Test Results and Report Templates

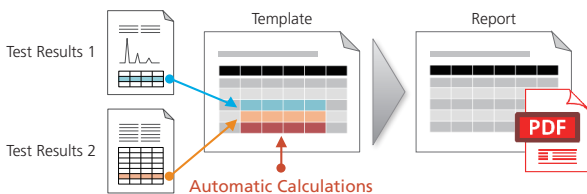

- ▶ Prevents manual transcription errors and tampering
- ▶ Shortens the processing time
- ▶ Prevents unauthorized alterations and deletions of reports


3 Automatic Transfer of Resulting Information to Upper Systems

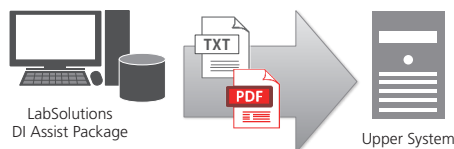

- ▶ Prevents manual transcription errors and tampering
- ▶ Shortens the processing time
- ▶ Establishes the format for output text files

Sample System Configuration

The LabSolutions DI Assist Package can be used in combination with a Shimadzu instrument. The PC is the same as the one that controls the Shimadzu instrument, so there is no need to prepare a new PC.

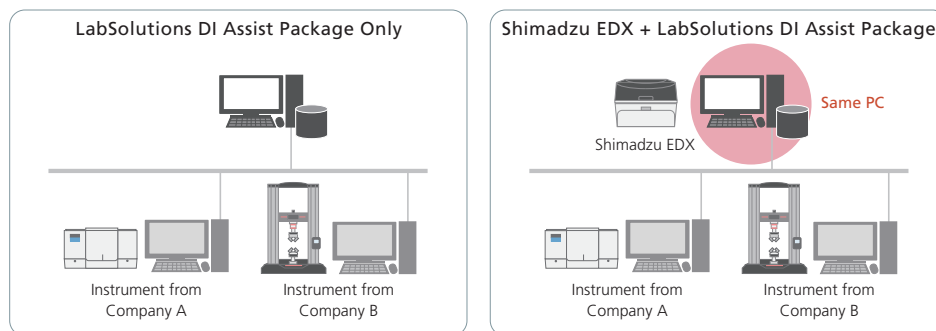


Note 1) Only one Shimadzu instrument can be connected to the PC for controlling Shimadzu instruments.

Note 2) Some instruments, including chromatography instruments, are not supported.

Expansion to LabSolutions™ CS

If there are a number of connected instruments and Shimadzu instruments, expansion to the LabSolutions CS network system is recommended. All information is managed on a server, heightening reliability and efficiency.

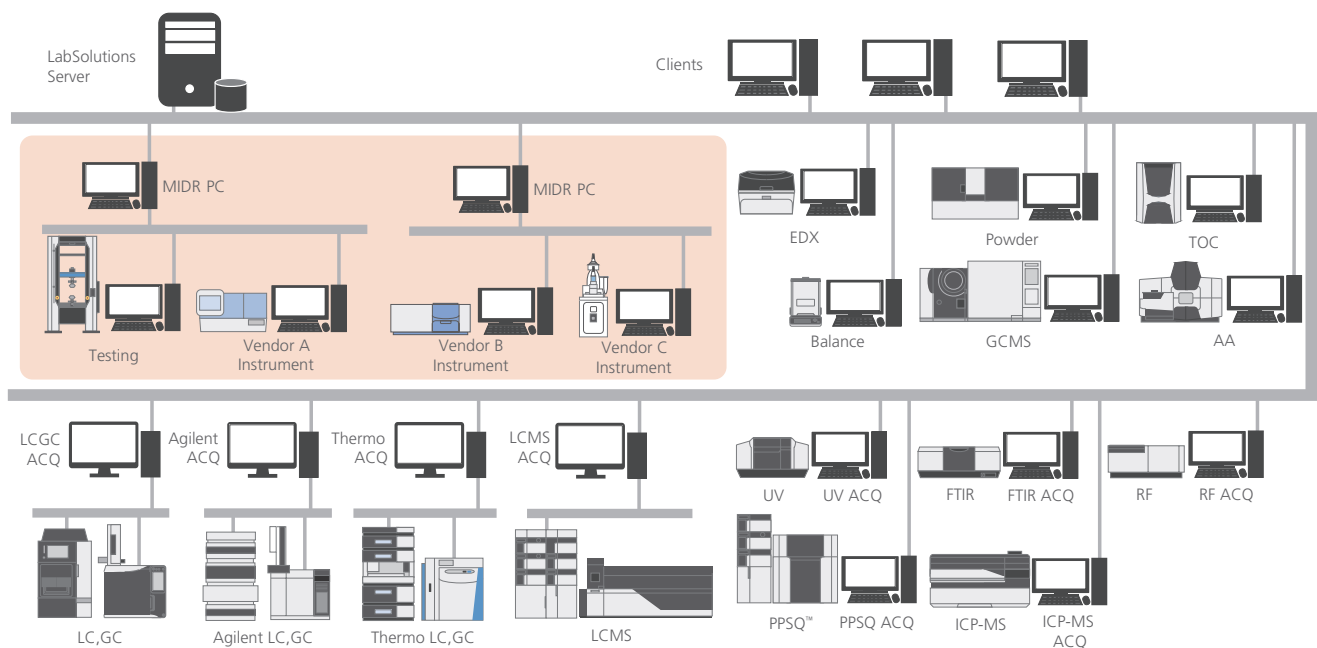

LabSolutions DI Assist Package, LabSolutions and PPSQ are trademarks of Shimadzu Corporation.



For Research Use Only. Not for use in diagnostic procedures.

This publication may contain references to products that are not available in your country. Please contact us to check the availability of these products in your country.

Company names, products/service names and logos used in this publication are trademarks and trade names of Shimadzu Corporation, its subsidiaries or its affiliates, whether or not they are used with trademark symbol "TM" or "®".

Third-party trademarks and trade names may be used in this publication to refer to either the entities or their products/services, whether or not they are used with trademark symbol "TM" or "®".

Shimadzu disclaims any proprietary interest in trademarks and trade names other than its own.

The contents of this publication are provided to you "as is" without warranty of any kind, and are subject to change without notice. Shimadzu does not assume any responsibility or liability for any damage, whether direct or indirect, relating to the use of this publication.

Shimadzu Corporation
www.shimadzu.com/an/