


InertSearch™ for GC

InertCap® Applications

Volatile organic compounds

Data No. GA005-0000

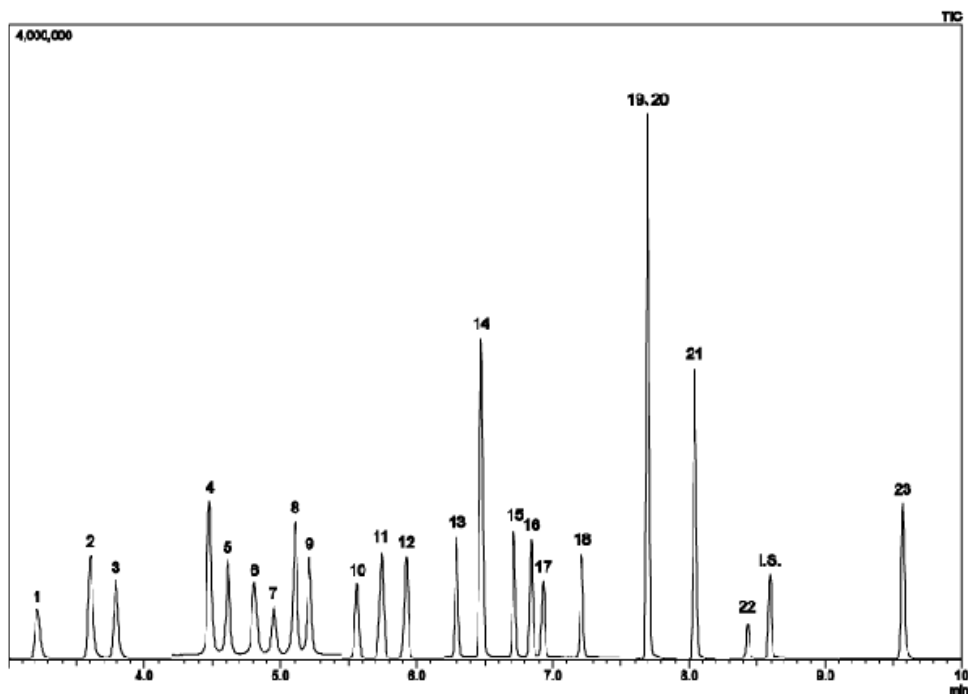

System : GC/MS P&T/GC

注) 高速分析用 HTP セットが必要になります。

Column : AQUATIC-H

Injection : P&T

Detection : MS

Solutes : 1. 1,1-Dichloroethene

2. Dichloromethane

3. t-1,2-Dichloroethene

4. cis-1,2-Dichloroethene

5. Chloroform

6. 1,1,1-Trichloroethane

7. Carbone Tetrachloride

8. 1,2-Dichloroethane

9. Benzene

10. Trichloroethene

11. 1,2-Dichloropropane

12. Bromodichloromethane

13. cis-1,3-Dichloro-1-propene

14. Toluene

15. t-1,3-Dichloro-1-propene

16. 1,1,2-Trichloroethane

17. Tetrachloroethene

18. Dibromochloromethane

19. m-Xylene

20. p-Xylene

21. o-Xylene

22. Bromoform

23. p-Dichlorobenzene

I.S. p-Bromofluorobenzene