


GC CAPILLARY COLUMN APPLICATIONS: ENVIRONMENTAL

Alcohols

COLUMN: BTR-CW, Carbowax/
BTR -Oxygen resistant
30M. x 0.53mm I.D. x
1.0µm film
Cat. No.: BTR-CW-30V-1.0F
Temperature: 50° (5°/min.) - 100°C
Injector: 220°C
Detector: 250°C , FID
Carrier Gas: 21.7 cm/sec., Helium

1. methanol
2. 2-propanol
3. ethanol
4. 2-methyl-2-butanol
5. 2-butanol
6. 1-propanol
7. 2-methyl-1-propanol
8. 3-pentanol
9. 2-pentanol
10. 1-butanol
11. 3-methyl-1-butanol
12. 1-pentanol

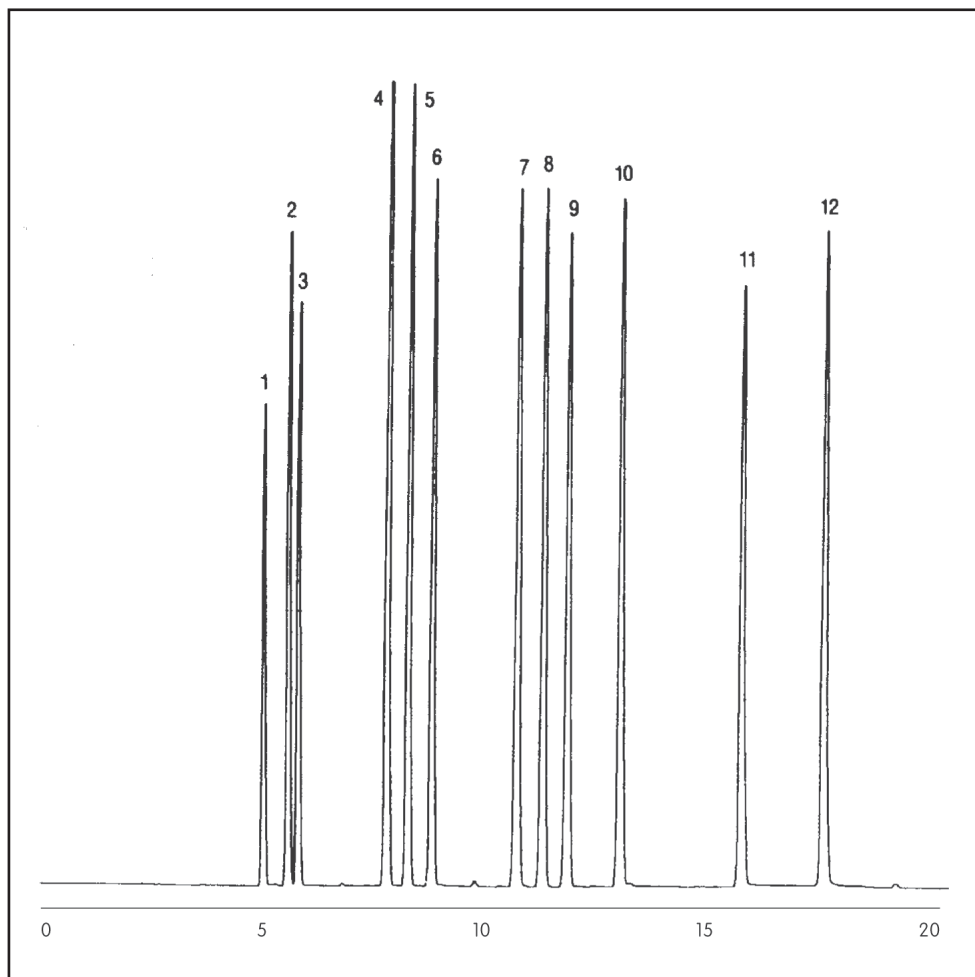

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com