


# GC CAPILLARY COLUMN APPLICATIONS: FAMES

## Fame Mix

COLUMN: 007-CW, Bonded  
Carbowax  
25M. x 0.25mm I.D. x  
0.25µm film  
Cat. No.: 007-CW-25-0.25F  
Temperature: 80° (8°/min.) - 240°C  
Injector: 230°C  
Detector: 280°C, FID  
Carrier Gas: 28 cm/sec., Helium

1. C-8:0
2. C-9:0
3. C-10:0
4. C-11:0
5. C-12:0
6. C-14:0
7. C-14:1
8. C-16:0
9. C-16:1
10. C-18:0
11. C-18:1
12. C-18:2
13. C-18:3
14. C-20:0
15. C-20:1
16. C-20:2
17. C-20:3

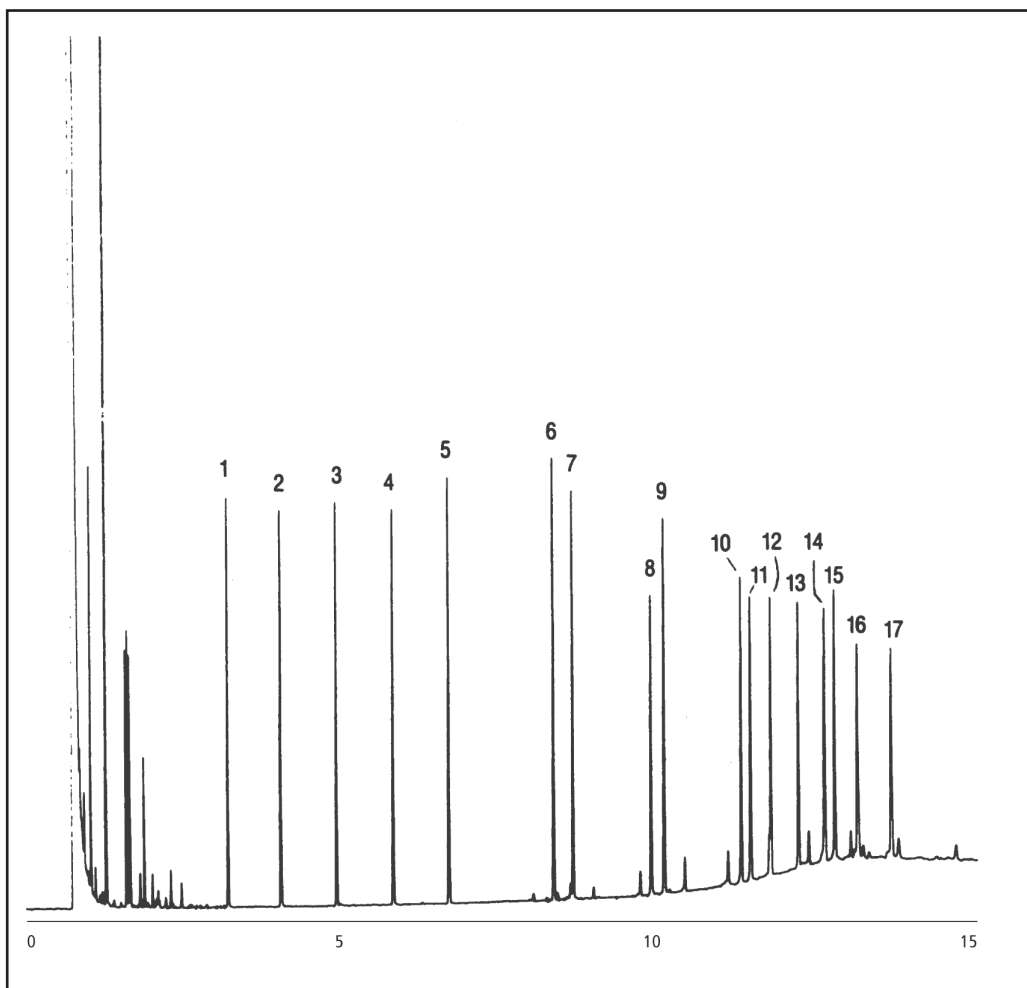

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391  
TOLLFREE: 800-275-7033 • E-MAIL: [tech@quadrexcorp.com](mailto:tech@quadrexcorp.com)